

external implant system



BIOHORIZONS[®]
SCIENCE • INNOVATION • SERVICE

Corporate Overview	1
External System Overview	2-3
Product Catalog	
Implants	4
Healing Abutments & Cover Screws	5
Surgical Kit & Components	6-7
Ancillary Surgical Instruments	8-9
Collateral Literature & Support Material	10-11
W&H Surgical Motors & Accessories	12-13
Abutment Overview	14
Impression Component Overview	15
Individual Abutments	16-19
Impression Components	20-21
Locator Abutment Components	22-23
Ball Abutment Components	24-25
Abutment for Screw Components	26-27
Prosthetic Instrumentation	28-29
Surgical Manual	
Introduction	30
Surgical Protocols	31
Implant Spacing	32
Surgical Kit Layout	33
Surgical Drills Overview	34
Osteotomy Initialization	35
Parallel Pins & Width Increasing Drills	36
Crestal Bone Drills & Bone Taps	37
Implant Packaging	38
Implant Pick-up	39
Abutment Removal	40
Cover Screws & Healing Abutments	41
Ancillary Instruments	42
Appendix	43
Icon Legend & References	44
Ordering & Warranty Information	45

BIOHORIZONS®

SCIENCE • INNOVATION • SERVICE

99.2%
average
implant
success
rate¹



BioHorizons is one of the fastest growing companies in the dental implant industry. Our comprehensive portfolio of dental implants and biologics products is evidence-based and scientifically-proven. From the launch of the External implant system (Maestro) in 1997, to the introduction of the Tapered Internal implant incorporating Laser-Lok® technology in 2007, BioHorizons has continued to provide clinicians unique products supported by university based research.

global
leader for
biologic
based
solutions



SCIENCE

BioHorizons uses science and innovation to create unique products with proven surgical and esthetic results.

INNOVATION

Our advanced implant technologies, biologic products and computer guided surgery software have made BioHorizons a leading dental implant company.

products
sold
in over
80
markets




SERVICE

BioHorizons understands the importance of providing excellent service. Our global network of professional representatives and our highly trained customer care support team are well-equipped to meet the needs of patients and clinicians.

external

Originally known as the Maestro, the External implant has performed masterfully in every published study since its launch in 1997. With success rates from 98.5% to 100% in studies of up to 1,400 implants in a wide variety of placement and loading conditions,¹ clinicians can confidently rely on the External to achieve consistent, reproducible results necessary to build a successful implantology practice.

- Transmits 10 times less destructive shear than conventional v-threads²
- Shown to achieve bone-to-implant contact levels of 80.6%³
- Higher reverse torque values than conventional v-thread implants⁴
- 4 platform/body diameters: 3.5mm, 4.0mm, 5.0mm, 6.0mm
- 5 lengths: 7mm, 9mm, 10.5mm, 12mm, 15mm

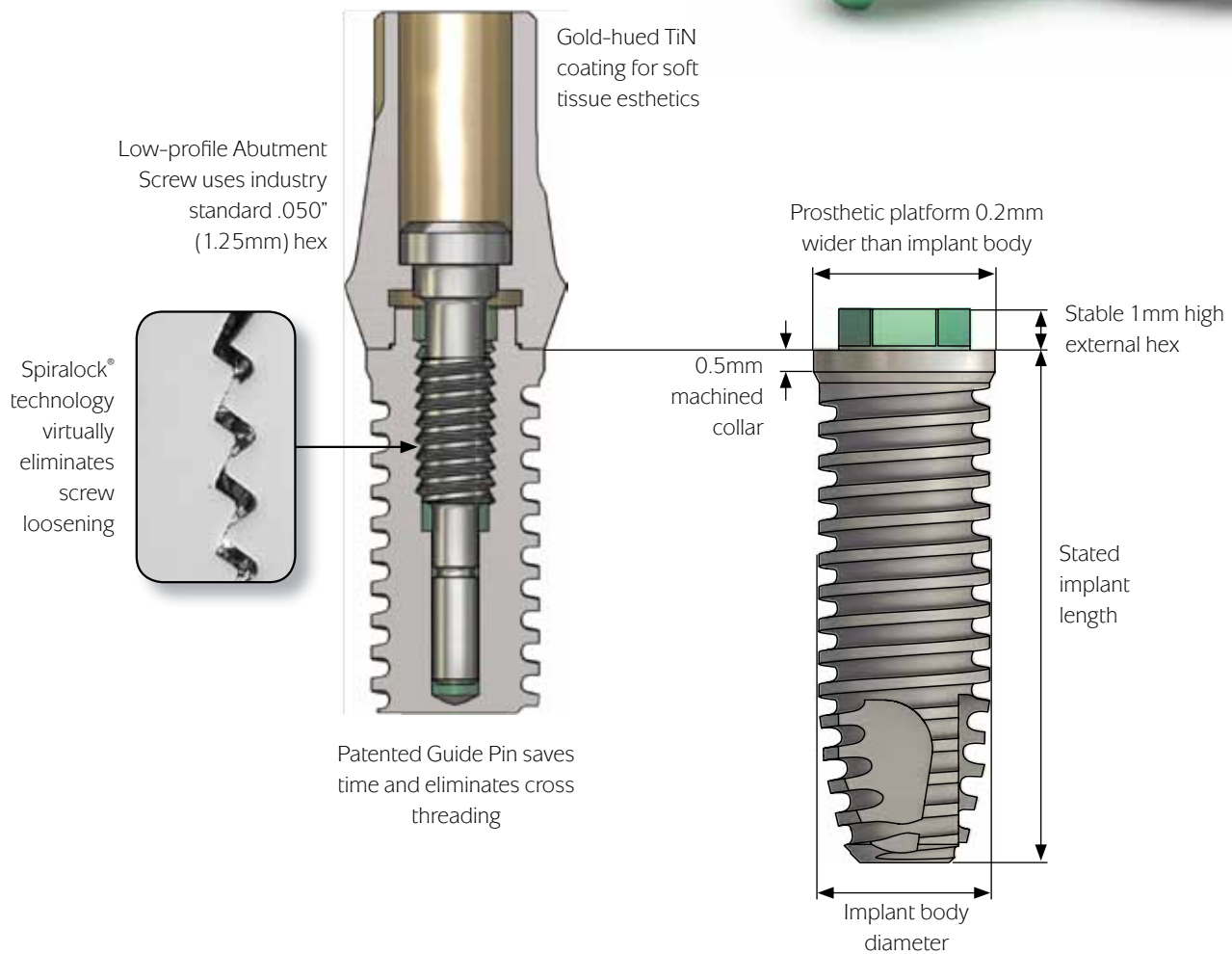

Crestal bone maintenance leads to success



A prospective, multicenter trial showed a 5-year cumulative success rate of 99.5% for 495 External implants in 151 patients. There were no significant differences in the results by center, implant type, bone density, area of the mouth or prosthesis type. Crestal bone loss from prosthesis delivery to one-year follow-up was extremely low (0.06mm), and actually turned into bone gain at later follow-up intervals.⁵


Radiographs of External implants in function for more than 6 years. Note that the crestal bone levels are maintained above the level of the first thread. Courtesy of Dr. Carl E. Misch.



Increased surface area for immediate loading success

Bone is strongest in compressive loading; weakest in shear. External's square-thread design imparts 10 times less destructive stresses at the implant / bone interface than conventional v-thread designs, while maximizing compressive load transfer and providing excellent primary stability. These attributes have become increasingly important with the emergence of early and immediate load protocols.^{6,7,8}

The square-thread design provides up to 154% greater surface area than conventional v-thread designs. A comparative study of thread designs (square-thread, v-thread and buttress thread) concluded that the square-thread design exhibited significantly higher reverse torque values and bone-to-implant contact than the other thread designs.³

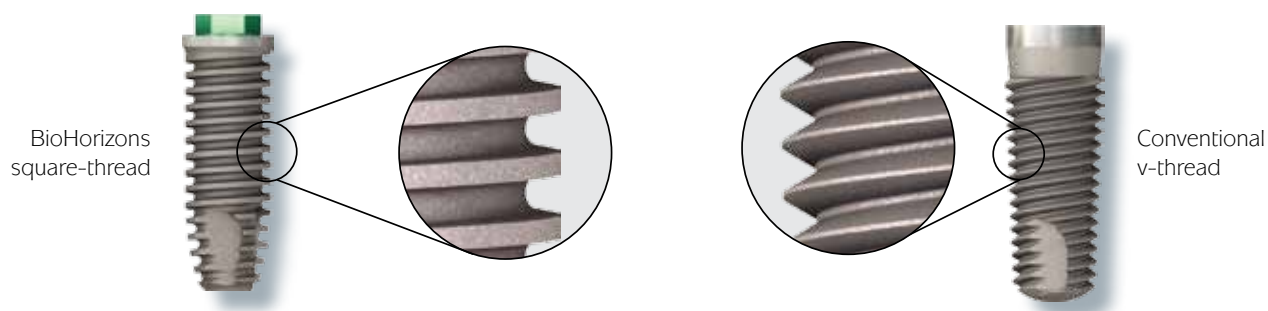

Biomechanical thread form | Backed by science

External implants feature the time-proven, biomechanical square thread which has been demonstrated to provide **increased bone-to-implant contact and yield higher reverse-torque values.**

	V-thread	Modified-square
Reverse Torque Value (Ncm)	15.58 ± 6.07	23.17 ± 9.68
% Bone-to-Implant Contact	65.46 ± 9.64	74.37 ± 8.63

Steigenga J, Al-Shammari K, Misch C, Nociti F and Wang H-L. *Effects of Implant Thread Geometry on Percentage of Osseointegration and Resistance to Reverse Torque in the Tibia of Rabbits.* J Periodontology 2004;75:1233-1241

Multiple studies have demonstrated the BioHorizons square thread to be extremely reliable in immediate loading protocols. One study, summarized below, followed 242 BioHorizons implants through immediate functional and non-functional loading.

BioHorizons Implants	N Implants	N Failures	% Implant Survival	% Prostheses Survival
Immediate Functional Loading	126	0	100%	100%
Immediate Non-Functional Loading	116	0	100%	100%

Degidi and Piattelli. *Immediate Functional and Non-Functional Loading of Dental Implants: A 2- to 60-Month Follow-Up Study of 646 Titanium Implants.* J Periodontology 2003;74:225-241.

EXTERNAL IMPLANTS

Resorbable Blast
Texturing (RBT)
Roughened
TiO₂ surface



Diameter x Length



Hydroxylapatite
Coating (HA)
High-density
crystalline coating

	3507D3 *	3.5mm x 7mm	NA	
	3509D3	3.5mm x 9mm	NA	
	35105D3	3.5mm x 10.5mm	NA	
	3512D3	3.5mm x 12mm	NA	
	3515D3	3.5mm x 15mm	NA	
	4007D3 *	4.0mm x 7mm	4007D4 *	
	4009D3	4.0mm x 9mm	4009D4	
	40105D3	4.0mm x 10.5mm	40105D4	
	4012D3	4.0mm x 12mm	4012D4	
	4015D3	4.0mm x 15mm	4015D4	
	5007D3 *	5.0mm x 7mm	5007D4 *	
	5009D3	5.0mm x 9mm	5009D4	
	50105D3	5.0mm x 10.5mm	50105D4	
	5012D3	5.0mm x 12mm	5012D4	
	5015D3	5.0mm x 15mm	5015D4	
	6007D3 *	6.0mm x 7mm	6007D4 *	
	6009D3	6.0mm x 9mm	6009D4	
	60105D3	6.0mm x 10.5mm	60105D4	
	6012D3	6.0mm x 12mm	6012D4	
	6015D3	6.0mm x 15mm	6015D4	

	Body Diameter	Platform Diameter	Apical Diameter	Hex flat-to-flat measure	Machined Collar Height	Minimum Ridge Width	Minimum Mesial/ Distal Space
3.5mm Implants	3.5mm	3.7mm	2.0mm	2.4mm	0.5mm	5.5mm	6.5mm
4.0mm Implants	4.0mm	4.2mm	2.1mm	2.7mm	0.5mm	6.0mm	7.0mm
5.0mm Implants	5.0mm	5.2mm	2.5mm	3.0mm	0.5mm	7.0mm	8.0mm
6.0mm Implants	6.0mm	6.2mm	3.5mm	3.0mm	0.5mm	8.0mm	9.0mm











Please refer to the External ~7mm Length Implants manual (ref. L0151) for detailed instructions for use of the prosthetic components for 7mm implants.

*Each 7mm implant comes packaged with the 3inOne Abutment, a short cover screw and a short Abutment Screw.

HEALING ABUTMENTS & COVER SCREWS

Healing Abutments

Used for tissue healing at second-stage uncover or as a transmucosal element for single-stage procedures. Color-coded by implant diameter. Hand-tighten with the .050" (1.25mm) Hex Driver. Titanium Alloy.

Flared Emergence		Implant Diameter x Height		Straight Emergence	
	212-303	3.5mm x 3mm	204-303		
	212-304	3.5mm x 4.5mm	204-304		
	212-306	3.5mm x 6mm	204-306		
	212-403	4.0mm x 3mm	204-403		
	212-404	4.0mm x 4.5mm	204-404		
	212-406	4.0mm x 6mm	204-406		
	212-503	5.0mm x 3mm	204-503		
	212-504	5.0mm x 4.5mm	204-504		
	212-506	5.0mm x 6mm	204-506		
	212-603	6.0mm x 3mm	204-603		
	212-604	6.0mm x 4.5mm	204-604		
	212-606	6.0mm x 6mm	204-606		

Cover Screws

Included with each implant system. Protects the implant's prosthetic platform during submerged surgical healing. Color-coded by implant diameter. Hand-tighten with the .050" (1.25mm) Hex Driver. Titanium Alloy.

	103-900	3.5mm Cover Screw
	104-900	4.0mm Cover Screw
	105-900	5.0mm Cover Screw
	106-900	6.0mm Cover Screw
	103-700	3.5mm Short Cover Screw
	104-700	4.0mm Short Cover Screw
	105-700	5.0mm Short Cover Screw
	106-700	6.0mm Short Cover Screw

products shown not to scale

External Surgical Kit



160-500

External Surgical Kit

Includes all instruments (except 122-12507) shown on pages 6-7.

160-200

Surgical Tray & Lid (without instruments)

Individual Surgical Kit Components



122-103
2.0mm Starter Drill



122-100
Drill Extender



144-100 Straight Parallel Pins (4 per kit)

144-200 20° Angled Parallel Pins (2 per kit)



122-12507 2.5 x 7mm Depth Drill (sold separately)

122-12509 2.5 x 9mm Depth Drill

122-125105 2.5 x 10.5mm Depth Drill

122-12512 2.5 x 12mm Depth Drill

122-12515 2.5 x 15mm Depth Drill

122-225 2.5mm Depth Drill (without Depth Stop)



122-230 3.0mm Width Increasing Drill

122-232 3.4mm Width Increasing Drill

122-237 3.9mm Width Increasing Drill

122-242 4.4mm Width Increasing Drill

122-247 4.9mm Width Increasing Drill

122-252 5.4mm Width Increasing Drill

SURGICAL KIT COMPONENTS


Individual Surgical Kit Components



- 123-000** 3.5mm Crestal Bone Drill
- 124-000** 4.0mm Crestal Bone Drill
- 125-000** 5.0mm Crestal Bone Drill
- 126-000** 6.0mm Crestal Bone Drill



- 123-300** 3.5mm Bone Tap
- 124-300** 4.0mm Bone Tap
- 125-300** 5.0mm Bone Tap
- 126-300** 6.0mm Bone Tap

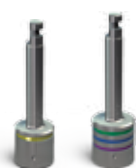

- 133-100** 3.5mm Handpiece Adapter
- 145-100** 4.0/5.0/6.0mm Handpiece Adapter



- 303-200** 3.5mm Ratchet Adapter*
- 345-200** 4.0/5.0/6.0mm Ratchet Adapter*



- 153-000** 3.5mm Implant-level Insertion Tool*
- 154-000** 4.0mm Implant-level Insertion Tool*
- 155-000** 5.0/6.0mm Implant-level Insertion Tool*



- 135-351**
- .050" (1.25mm) Hex Driver**



- 130-000**
- Ratchet**



- 300-400**
- Hand Wrench**



- 300-206**
- 4mm Square Drive Extender***
- Replaced 300-205 starting in June 2010.*
- Includes PEEK C-ring for durable retention in Ratchet. Cannot be used with bone taps.*



- 144-300**
- Implant Spacer / Depth Probe**

**instrument o-rings & c-rings wear out over time. If an instrument is no longer held securely by its associated driver, order a replacement ring through Customer Care.
products shown not to scale*

Extended Shank Drills *(not included in Surgical Kit)*

Extended Shank Drills are 8mm longer than our standard drills. They provide an intermediary increase in length between the standard drills used alone, or with the Drill Extender (ref. 122-100), which adds 16mm of length. The depth marks and cutting geometry are identical to our standard drills.



122-403	2.0mm Starter Drill, Extended Shank
122-42507	2.5 x 7mm Depth Drill, Extended Shank
122-42509	2.5 x 9mm Depth Drill, Extended Shank
122-425105	2.5 x 10.5mm Depth Drill, Extended Shank
122-42512	2.5 x 12mm Depth Drill, Extended Shank
122-42515	2.5 x 15mm Depth Drill, Extended Shank
122-425	2.5mm Depth Drill, Extended Shank



122-430	3.0mm Width Increasing Drill, Extended Shank
122-432	3.4mm Width Increasing Drill, Extended Shank
122-437	3.9mm Width Increasing Drill, Extended Shank
122-442	4.4mm Width Increasing Drill, Extended Shank
122-447	4.9mm Width Increasing Drill, Extended Shank
122-452	5.4mm Width Increasing Drill, Extended Shank

Ratchet & Hand Wrench Extender



300-205	Ratchet & Hand Wrench Extender
----------------	---

Can be used with Bone Taps (page 7).

Surgical Driver



150-000	Surgical Driver
----------------	------------------------

Used to drive implants into the osteotomy, particularly in the anterior region. The driver holds the Abutment-level Driver, Ratchet which interfaces with the *3inOne* Abutment. Also interfaces with the .050" (1.25mm) Hex Drivers as well as Bone Taps and the Implant-level Drivers, Ratchet.

ANCILLARY SURGICAL INSTRUMENTS

Ancillary Instruments *(not included in Surgical Kit)*



122-110 2.0mm Lindemann Bone Cutter

Side-cutting drill used to correct eccentric osteotomy preparations.



PYTP 3.5mm Tissue Punch

122-201 4.0mm Tissue Punch

122-202 5.0mm Tissue Punch

122-203 6.0mm Tissue Punch

BioHorizons Tissue Punches are used in a latch-type handpiece to remove the soft tissue from the crest of the ridge prior to osteotomy preparation in a flapless surgical procedure.



122-113 External Bone Profiler, 3.5mm

122-114 External Bone Profiler, 4.0mm

122-115 External Bone Profiler, 5.0mm

122-116 External Bone Profiler, 6.0mm

Used at implant uncoverly to contour crestal bone to accommodate abutments when the implant is subcrestal. For use in latch-type reduction handpieces. The Profiler's internal geometry matches the geometry of the surgical Cover Screw which aligns the Profiler for precise removal of tissue surrounding the platform.



134-350 .050" (1.25mm) Handpiece Hex Driver

134-450 .050" (1.25mm) Handpiece Hex Driver, Long

For installation and removal of Cover Screws, Healing Abutments and Abutment Screws. Use with latch-type contra-angle handpieces: WS-75 E/KM or similar. 134-450 is 5mm longer than the 134-350.

Adjustable Torque Wrenches



ATW ITL Precise Adjustable Torque Wrench

Ratchet design places both implants and abutments with 9 distinct torque settings (15, 20, 25, 30, 35, 40, 45, 50 and 60 Ncm). A simple twist of the handle locks in precision-engineered torque values and guarantees accuracy and repeatability.




AGYR-15500 Anthogyr Torq Control® Torque Wrench

Ergonomic rotating handle design is the ideal solution for access to screws placed in the posterior. The 7 predetermined torque values (10, 15, 20, 25, 30, 32 and 35 Ncm) make it a tool of extreme precision.

products shown not to scale


Related Materials

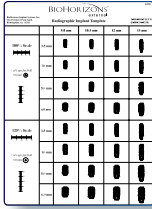

L0250

External Radiographic Implant Template (overlay)

Designed to aid the clinician in the preoperative determination of options for implant length and diameter. The clear overlay template shows all sizes of BioHorizons External implants in 100% and 125% scale. Packaged with each new surgical kit, but also available upon request.

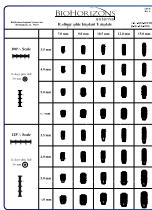

L0136

External Radiographic Implant Template w/7mm

Designed to aid the clinician in the preoperative determination of options for implant length and diameter. The clear overlay template shows all sizes of BioHorizons External implants in 100% and 125% scale.



L0151

External - 7mm Length Implants

Guides the clinician and laboratory through the nuances of placing and restoring BioHorizons External 7mm length implants. Includes ordering information for implants, 2.0 x 7mm Depth Drills and Radiographic Templates with 7mm length implants.


NOTE: 7mm length implants are not available in all countries; please call for availability in your market.



L0130

External Surgical Kit Wall Chart

11" x 17" (28cm x 43cm) wall chart details the location of instruments in the External Surgical Kit (ref. 160-500), as well as the reference numbers to assist re-ordering. Packaged with each new surgical kit, but also available upon request.

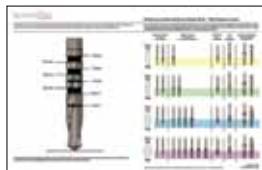

L0131

External Depth Mark and Drill Sequence Chart

11" x 17" (28cm x 43cm) wall chart details the Drill Depth Marks and recommended sequence of use for the instruments found in the External Surgical Kit (ref. 160-500). Packaged with each new surgical kit, but also available on request.



ML0207

External Impression Technique Guide

Guides the clinician and laboratory through impression making techniques and procedures with BioHorizons External impression components.

Patient Education Materials

ML0103 **Patient Education – Tooth replacement with dental implants**

ML0114 **Patient Education – Stabilizing dentures using dental implants**

ML0129 **Patient Education – Rebuilding and maintaining bone**

MLD101 **Patient Education – Soft tissue grafting with AlloDerm**



ML0131 **Patient Education – Dental Implants - the tooth replacement solution**

This high-quality flipbook helps the implant candidate understand the rationale and the advantages of implant therapy compared to traditional treatment methods. 9" x 6" (23cm x 16cm).



M4XIM **External 4x Patient Education Model**

This 4-times scale model includes the Implant, Abutment, Abutment Screw, Ball-top Screw, Healing Abutment, Cover Screw, Implant Driver and Hex Driver. Excellent for demonstrating implant components to the patient. Call for availability.



MPSA **Five Implant Acrylic Model**

This life-sized model allows the clinician to illustrate implant placement of BioHorizons implants. Clear acrylic allows the implants to be viewed in relation to adjacent teeth. Call for availability.



EP-MSLA **Locator Patient Education Model**

This life-sized model allows the clinician to illustrate a denture on an edentulous mandible supported by 2 implants with Locator abutments. Clear acrylic allows the implants to be viewed in relation to the denture. Call for availability.



SAMPLE 4012D3 **Sample External Implant**

Sample implants are good for implant demonstrations and sawbones labs. These sample implants are fully anodized and are not for human use. Available in 4.0 x 12mm, D3 threadform.





BioHorizons proudly distributes W&H implant motors, handpieces and accessories. Additional W&H products and re-order items are available. For more information, contact your BioHorizons representative or visit the online catalog (www.biohorizons.com).

W&H ElcoMED SA-310 Professional Kits

Professional Kits include: ElcoMED SA-310 console, motor with 1.8m attached cable, premium foot pedal, (9) complete disposable irrigation tubes, USB stick, handpiece, bur testing gauge, spray cap and service oil (ships separately).



WH-310L

ElcoMED SA-310 Professional Kit with LED

Includes LED handpiece (WH-10207530).

WH-310

ElcoMED SA-310 Professional Kit

Includes handpiece (WS-75).

W&H ImplantMED SI-915 Starter Kits

Starter Kits include: ImplantMED SI-915 console, motor with 1.8m attached cable, standard foot pedal, (9) complete disposable irrigation tubes, handpiece, bur testing gauge, spray cap and service oil (ships separately).



WH-915L

ImplantMED 915 Starter Kit with LED

Includes mono block LED handpiece (WH-10207560).

WH-915

ImplantMED SI-915 Starter Kit

Includes mono block handpiece (WI-75MB).

W&H ElcoMED SA-310 and ImplantMED SI-915 re-order items

Ref. Code	Description
WH-ITC915	Disposable irrigation tubing, 2.2m ImplantMED SI-915 and ElcoMED SA-310 (box of 6)
WH-IC	Irrigation spray clip for external and internal irrigation (set of 3)
WH-MD400	MD-400 service-oil F1
WH-OSC	Oil spray cap for WH-MD400

Surgical Handpieces

WH-10207530

WS-75 E/KM LED G Surgical 20:1 Contra-Angle


LED, 20:1, hexagon chucking system, press-button chuck system, for surgical burs and cutters with a 2.35mm contra-angle shank, internal and external cooling system. Dismantle for easy cleaning.



WS-75

WS-75 E/KM Surgical 20:1 Contra-Angle

20:1, hexagon chucking system, press-button chuck system, for surgical burs and cutters with a 2.35mm contra-angle shank, internal and external cooling system. Dismantle for easy cleaning.

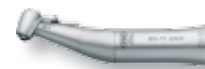

WH-10207560

WI-75 E/KM LED G Surgical 20:1 Contra-Angle, Mono Block

LED, 20:1, hexagon chucking system, press-button chuck system, for surgical burs and cutters with a 2.35mm contra-angle shank, internal and external cooling system.



WI-75MB

WI-75 E/KM Surgical 20:1 Contra-Angle, Mono Block

20:1, hexagon chucking system, press-button chuck system, for surgical burs and cutters with a 2.35mm contra-angle shank, internal and external cooling system.



WS-92

WS-92 E/3 Surgical 1:2.7 Contra-Angle

1:2.7, press-button chuck system with triple spray, for surgical burs and cutters with a 1.6mm friction grip shank. Dismantle for easy cleaning.



S-12

S-12 Surgical 1:2 Angled

1:2, angled, extra slim tip for increased vision, lever chuck system, for surgical burs and cutters 2.35mm >70mm extra long shaft. Dismantle for easy cleaning.



S-11

S-11 Surgical 1:1 Straight

1:1, straight, lever chuck system, for surgical burs and cutters 2.35mm >45mm shaft. Dismantle for easy cleaning.



Prosthodontic Screw Driver

WH-16934000

IA-400 Prosthodontic Screw Driver

Handpiece drive with key pad, 80:1 contra-angle handpiece with hexagon chucking system, charging station, rechargeable Li-ion battery, and power cable.



Bur Testing Gauge

WH-BTG


Bur Testing Gauge

Used to verify the condition of latch-type burs. Burs in proper condition will fit into larger diameter hole, but will not fit into the smaller hole (marked red). Burs that fail either of these tests are unfit for use, and may cause damage to the handpiece.




IMPLANT-LEVEL FIXED ABUTMENT OVERVIEW

Cement-retained

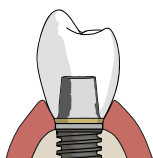


Abutments for Cement-retention
pages 16-17

Custom fabricated abutments

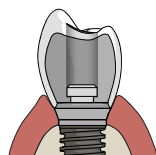

Custom Castable Abutments
page 18


Provisional Restorations

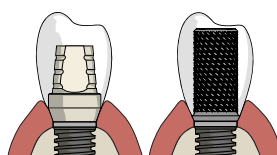

Temporary Abutments
page 19



REMOVABLE AND HYBRID ABUTMENT OVERVIEW

Implant-level Cast Bar or Hybrid (Fixed-detachable)

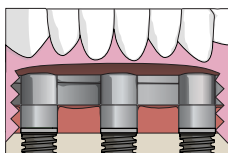


Custom Castable Abutments
page 18


Abutment for Screw level Cast Bar or Hybrid (Fixed-detachable)

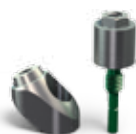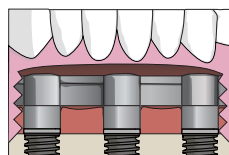

Abutments for Screw
pages 26-27

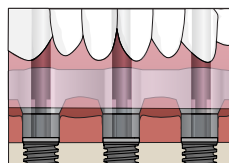


Abutment-retained Overdenture, Tissue-supported

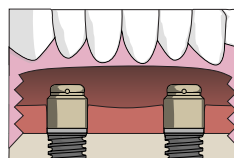


Locator Abutments
& components
pages 22-23

Ball Abutments
& components
pages 24-25

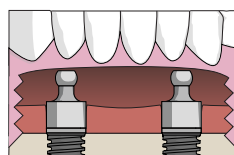


IMPRESSION COMPONENTS FOR IMPLANT-LEVEL RESTORATIONS

Implant-level Closed-tray Transfer

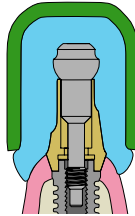


Indirect Transfer Copings
page 19

Implant-level Open-tray Pick-up

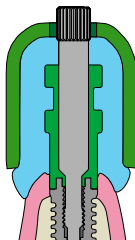

Direct Pick-up Copings
pages 20


Direct Abutment Impression (*Crown & Bridge-type*)

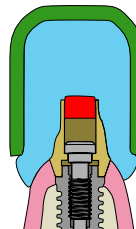

Abutments for cement-retention
pages 16-17


IMPRESSION COMPONENTS FOR ABUTMENT-LEVEL RESTORATIONS

Abutment for Screw-level Closed-tray Transfer

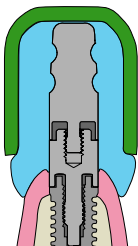


Abutment for Screw-level Indirect
Transfer Copings
page 26-27

Abutment for Screw-level Open-tray Pick-up

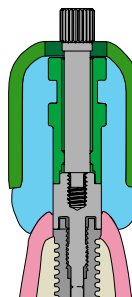


Abutment for Screw-level Direct
Pick-up Copings
page 26-27

Locator Abutment Closed-tray Pick-up

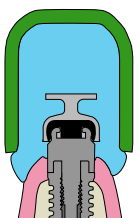


Locator components
pages 22-23

Ball Abutment Closed-tray Transfer

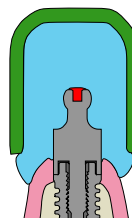

Ball Abutment components
pages 24-25


Ceramic Abutments for Esthetic Zone Restorations

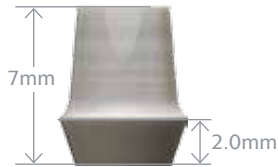

- 223-315** 3.5mm Ceramic Abutment
- 224-315** 4.0mm Ceramic Abutment
- 225-315** 5.0mm Ceramic Abutment
- 226-315** 6.0mm Ceramic Abutment

Used to fabricate cement-retained, single- or multiple-unit prostheses. Packaged with Abutment Screw (130-400 or 140-400 used only with External Ceramic Abutments). Yttria-stabilized Zirconia. Final torque: 30Ncm.


3inOne Abutments for Cement-retained Restorations

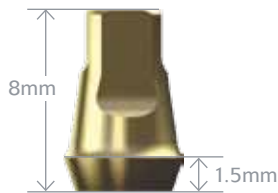

- 203-311** 3.5mm 3inOne Abutment
- 204-411** 4.0mm 3inOne Abutment
- 205-511** 5.0mm 3inOne Abutment
- 206-611** 6.0mm 3inOne Abutment

Used to fabricate cement-retained, single- or multiple-unit prostheses. Also used with Ball-top Screw as a closed-tray, hexed-timed impression coping. Packaged with Abutment Screw (130-300 or 140-300). Titanium Alloy. TiN coated. Final torque: 30Ncm.


Angled Abutments for Cement-retained Restorations

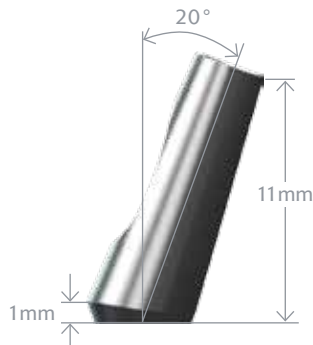

- 213-400** 3.5mm Angled Abutment
- 214-400** 4.0mm Angled Abutment
- 215-400** 5.0mm Angled Abutment

Used to fabricate cement-retained, single- or multiple-unit prostheses. Packaged with Abutment Screw (130-300 or 140-300). Titanium Alloy. Final torque: 30Ncm.

Angled Esthetic Abutments for Cement-retained Restorations

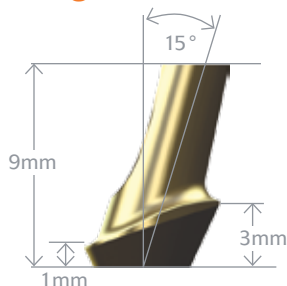

- 223-313** 3.5mm Angled Esthetic Abutment
- 224-313** 4.0mm Angled Esthetic Abutment
- 225-313** 5.0mm Angled Esthetic Abutment


Used to fabricate cement-retained, single- or multiple-unit prostheses. Packaged with Abutment Screw (130-300 or 140-300). Titanium Alloy. TiN coated. Final torque: 30Ncm.

INDIVIDUAL ABUTMENTS

Straight Abutments for Cement-retained Restorations

- 233-410** 3.5mm Straight Abutment
- 234-410** 4.0mm Straight Abutment
- 235-410** 5.0mm Straight Abutment


Used to fabricate cement-retained, single- or multiple-unit prostheses. Packaged with Abutment Screw (130-300 or 140-300). Titanium Alloy. Final torque: 30Ncm.

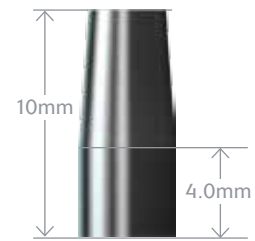

Straight Esthetic Abutments for Cement-retained Restorations

- 223-314** 3.5mm Straight Esthetic Abutment
- 224-314** 4.0mm Straight Esthetic Abutment
- 225-314** 5.0mm Straight Esthetic Abutment


Used to fabricate cement-retained, single- or multiple-unit prostheses. Packaged with Abutment Screw (130-300 or 140-300). Titanium Alloy. TiN coated. Final torque: 30Ncm.

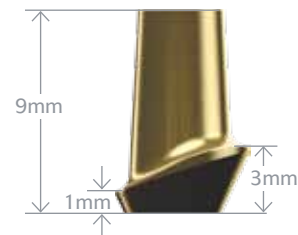

Profile Abutments for Cement-retained Restorations

- 223-312** 3.5mm Profile Abutment
- 224-312** 4.0mm Profile Abutment
- 225-312** 5.0mm Profile Abutment

Used to fabricate cement-retained, single- or multiple-unit prostheses. Packaged with Abutment Screw (130-300 or 140-300). Titanium Alloy. Final torque: 30Ncm.

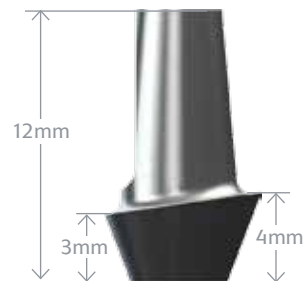


Abutments for Cement-retained Restorations - Non-hexed

- 223-208** 3.5mm Abutment for Cement, Non-hexed
- 224-208** 4.0mm Abutment for Cement, Non-hexed
- 225-208** 5.0mm Abutment for Cement, Non-hexed


Used to fabricate cement-retained, multiple-unit prostheses. One-piece, non-hex engaging abutment. Titanium Alloy. Final torque: 30Ncm.

Note: Must be prepped chairside followed by a conventional crown and bridge impression.

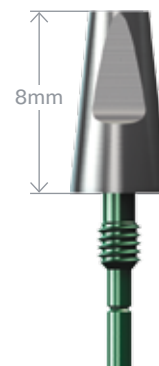

products shown not to scale

Custom Castable (UCLA) Abutments - Hexed

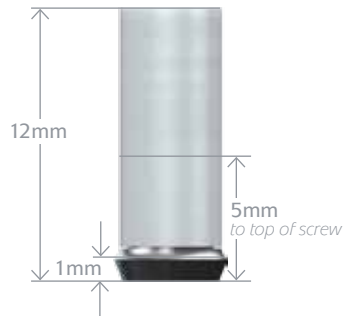


233-227 3.5mm Gold/Plastic Custom Castable Abutment, Hexed

234-227 4.0mm Gold/Plastic Custom Castable Abutment, Hexed

235-227 5.0mm Gold/Plastic Custom Castable Abutment, Hexed

236-227 6.0mm Gold/Plastic Custom Castable Abutment, Hexed

Used for single-unit screw-retained or custom abutment cement-retained restorations. Packaged with the Abutment Screw (130-300 or 140-300). Gold Alloy base with natural acetyl (Delrin® or Pomalux®) sleeve. Final torque: 30Ncm.

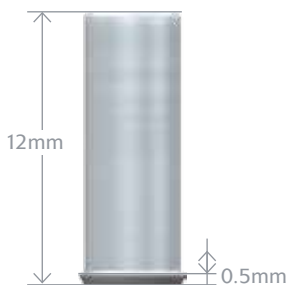

233-217 3.5mm Plastic Custom Castable Abutment, Hexed


234-217 4.0mm Plastic Custom Castable Abutment, Hexed

235-217 5.0mm Plastic Custom Castable Abutment, Hexed

236-217 6.0mm Plastic Custom Castable Abutment, Hexed

Used for multiple-unit, screw-retained restorations. Packaged with the Abutment Screw (PXAS). Natural acetyl (Delrin® or Pomalux®) sleeve. Final torque: 30Ncm.

Custom Castable (UCLA) Abutments - Non-hexed

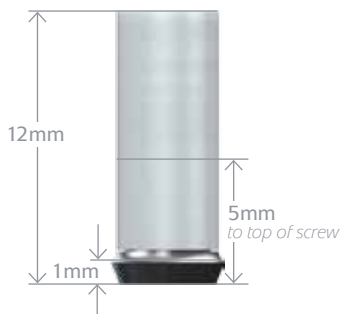


233-228 3.5mm Gold/Plastic Custom Castable Abutment, Non-hexed

234-228 4.0mm Gold/Plastic Custom Castable Abutment, Non-hexed

235-228 5.0mm Gold/Plastic Custom Castable Abutment, Non-hexed

236-228 6.0mm Gold/Plastic Custom Castable Abutment, Non-hexed

Used for single-unit screw-retained or custom abutment cement-retained restorations. Packaged with the Abutment Screw (130-300 or 140-300). Gold Alloy base with natural acetyl (Delrin® or Pomalux®) sleeve. Final torque: 30Ncm.

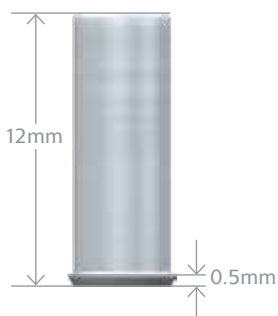

233-210 3.5mm Plastic Custom Castable Abutment, Non-hexed

234-210 4.0mm Plastic Custom Castable Abutment, Non-hexed

235-210 5.0mm Plastic Custom Castable Abutment, Non-hexed

236-210 6.0mm Plastic Custom Castable Abutment, Non-hexed


Used for multiple-unit, screw-retained restorations. Packaged with the Abutment Screw (PXAS). Natural acetyl (Delrin® or Pomalux®) sleeve. Final torque: 30Ncm.

INDIVIDUAL ABUTMENTS

PEEK® Temporary Abutments

- 273-000** 3.5mm PEEK Temporary Abutment
- 274-000** 4.0mm PEEK Temporary Abutment
- 275-000** 5.0mm PEEK Temporary Abutment
- 276-000** 6.0mm PEEK Temporary Abutment


For fabrication of cement- or screw-retained provisional restorations (up to 30 days). A Direct Coping Screw (243-320 or 244-320, purchased separately) may be used to maintain screw access hole during fabrication of screw-retained provisional prostheses. Packaged with Abutment Screw (130-300 or 140-300). PEEK (PolyEtherEtherKetone). Final torque: 30Ncm.

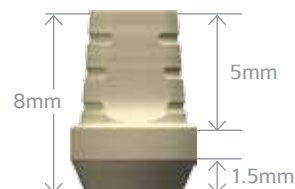

Titanium Temporary Abutments - Hexed

- 233-215** 3.5mm Titanium Temporary Abutment, Hexed
- 234-215** 4.0mm Titanium Temporary Abutment, Hexed
- 235-215** 5.0mm Titanium Temporary Abutment, Hexed


Used for single-unit screw-retained, long term temporary restorations (>30 days). Packaged with the Abutment Screw (130-300 or 140-300, hand-tighten). Titanium Alloy. Final torque: 30Ncm.

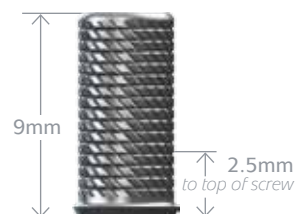

Titanium Temporary Abutments - Non-hexed

- 233-214** 3.5mm Titanium Temporary Abutment, Non-hexed
- 234-214** 4.0mm Titanium Temporary Abutment, Non-hexed
- 235-214** 5.0mm Titanium Temporary Abutment, Non-hexed

Used for multiple-unit, screw-retained, long term temporary restorations (>30 days). Packaged with the Abutment Screw (130-300 or 140-300, hand-tighten). Titanium Alloy. Final torque: 30Ncm.

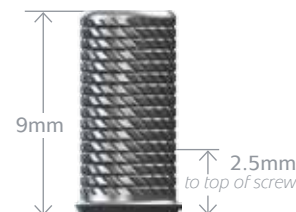

Ball-top Screw for Indirect Transfer


- 213-235** 3.5mm Ball-top Screw
- 214-235** 4.0/5.0/6.0mm Ball-top Screw

213-235 Fits 3.5mm prosthetic platforms. 214-235 fits 4/5/6.0mm prosthetic platforms. Used with the *3inOne* Abutment to form an impression coping for closed-tray, hexed-timed transfers. Hand-tighten. Titanium Alloy.



products shown not to scale

Direct Pick-up Copings (Open Tray) - Hexed

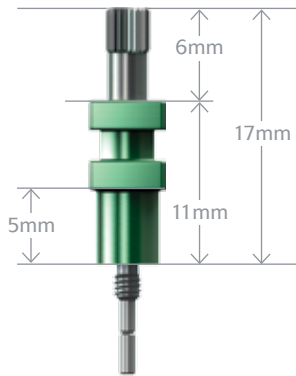

- 243-331** 3.5mm Direct Pick-up Coping, Hexed
- 244-331** 4.0mm Direct Pick-up Coping, Hexed
- 245-331** 5.0mm Direct Pick-up Coping, Hexed
- 246-331** 6.0mm Direct Pick-up Coping, Hexed

Used to make an open-tray, implant-level, hexed-timed impression. Packaged with the extended Direct Coping Screw (243-320 or 244-320). Titanium Alloy. Hand-tighten.

Direct Pick-up Copings (Open Tray) - Non-hexed



- 243-321** 3.5mm Direct Pick-up Coping, Non-hexed
- 244-321** 4.0mm Direct Pick-up Coping, Non-hexed
- 245-321** 5.0mm Direct Pick-up Coping, Non-hexed
- 246-321** 6.0mm Direct Pick-up Coping, Non-hexed

Used to make an open-tray, implant-level, non-hexed impression. Packaged with the extended Direct Coping Screw (243-320 or 244-320). Titanium Alloy. Hand-tighten.

Implant Analogs



- 293-000** 3.5mm Implant Analog
- 294-000** 4.0mm Implant Analog
- 295-000** 5.0mm Implant Analog
- 296-000** 6.0mm Implant Analog

Represents the implant in a laboratory-fabricated, implant-level stone model. Color-coded by implant body diameter. Titanium Alloy.

Direct Coping Screws



243-320 3.5mm Direct Pick-up Screw

244-320 4.0/5.0/6.0mm Direct Pick-up Screw

243-320 fits 3.5mm prosthetic platforms. 244-320 fits 4/5/6.0mm prosthetic platforms. Packaged with all Direct Pick-up Copings. May also be used in place of the Abutment Screw when extra length is needed, or to maintain the screw access hole during fabrication of screw-retained provisional prostheses. Up to 7mm can be removed from the screw without losing the hex engagement. Utilizes the .050" (1.25mm) Hex Driver. Hand-tighten or torque to 30 Ncm depending on application. Titanium Alloy.

Abutment Screws



130-300 3.5mm Abutment Screw

140-300 4.0/5.0/6.0mm Abutment Screw

Low profile screw head. Packaged with all two-piece abutments. Utilizes the .050" (1.25mm) Hex Driver. Final torque: 30 Ncm. Titanium Alloy.

Abutment Screws for 7mm length implants





130-050 3.5mm Short Abutment Screw, for 7mm

140-050 4/5/6mm Short Abutment Screw, for 7mm

Specifically designed for and packaged with each 7mm length External Implant system. May be used with 7mm length implants "as is"; or with "off-the-shelf" components *in lieu* of the regular Abutment Screw. Utilizes .050" (1.25mm) Hex Driver. Titanium Alloy.

Locator® Abutments

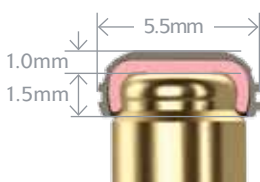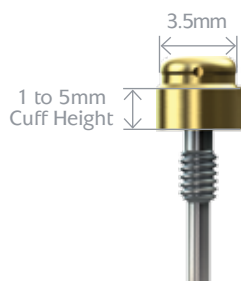

LA-8531	3.5mm Locator Abutment, 1.25mm Cuff Height
LA-8532	3.5mm Locator Abutment, 2.0mm Cuff Height
LA-8533	3.5mm Locator Abutment, 3.0mm Cuff Height
LA-8534	3.5mm Locator Abutment, 4.0mm Cuff Height
LA-8535	3.5mm Locator Abutment, 5.0mm Cuff Height

LA-8541	4.0mm Locator Abutment, 1.0mm Cuff Height
LA-8542	4.0mm Locator Abutment, 2.0mm Cuff Height
LA-8543	4.0mm Locator Abutment, 3.0mm Cuff Height
LA-8544	4.0mm Locator Abutment, 4.0mm Cuff Height
LA-8545	4.0mm Locator Abutment, 5.0mm Cuff Height

LA-8551	5.0mm Locator Abutment, 1.5mm Cuff Height
LA-8552	5.0mm Locator Abutment, 2.5mm Cuff Height
LA-8553	5.0mm Locator Abutment, 3.5mm Cuff Height
LA-8554	5.0mm Locator Abutment, 4.5mm Cuff Height
LA-8555	5.0mm Locator Abutment, 5.5mm Cuff Height

The Locator Implant Attachment System is designed for use with overdentures or partial dentures retained in whole or in part by dental implants in the mandible or maxilla. Order by Cuff Height to match the height of the gingival tissue. The abutment will extend above the tissue by 1.5mm to allow the Locator Male to seat completely. Order one Locator Male Processing Set for each Locator Abutment (sold in packs of 2 or 10). Locator Abutments are made from Titanium Alloy.

The Male Processing Package provides 3 choices of retention (see opposite page). The Replacement Males (clear, pink and blue) are used to restore implants with up to 10° of divergence (20° between implants). The Extended Range Replacement Males (green and red) accommodate divergences from 10° and 20° (40° between implants), and may be purchased separately.

Locator Components

LCT



Core Tool

Multi-purpose tool serves as hand driver for seating Locator Abutments onto the implants, seating tool for nylon male inserts and insert removal tool.

LMPP-2



Locator Male Processing Package (2 pack)

Includes: (2) Denture Caps with (2) Black Processing Males; (2) White Block-out Spacers; (2) Clear, (2) Pink and (2) Blue Nylon Males.

LMPP-10

Locator Male Processing Package (10 pack)

Includes: (10) Denture Caps with (10) Black Processing Males; (10) White Block-out Spacers; (10) Clear, (10) Pink and (10) Blue Nylon Males.

LOCATOR ABUTMENT COMPONENTS

Locator Impression Coping

LIC **Locator Impression Coping (4 pack)**



Locator Female Analog

LFA-4MM **Locator Female Analog 4mm (4 pack)**
Used for 3.5/4.5 platform

LFA-5MM **Locator Female Analog 5mm (4 pack)**
Used for 5.7 platform



Locator Male Replacements

LRM-C **Locator Replacement Male (clear) (4 pack)**
Retention: 5lb / 2268g

LRM-P **Locator Light Retention Replacement Male (pink) (4 pack)**
Retention: 3lb / 1361g

LRM-B **Locator Extra Light Retention Replacement Male (blue) (4 pack)**
Retention: 1.5lb / 680g

LRM-G **Locator Extended Range Replacement Male (green) (4 pack)**
Retention: 3-4lbs / 1361-1814g

LRM-R **Locator Extended Range Extra Light Retention Rep. Male (red) (4 pack)**
Retention: 1.5lbs / 680g

LBPRM **Locator Black Processing Replacement Male (4 pack)**



Locator Parallel Post

LPP **Locator Parallel Post (4 pack)**



Locator Measurement Guide

LAMG **Locator Angle Measurement Guide**



Locator Square Drive Tool

LSDT-15MM **Locator Square Drive Tool (15mm)**
Used with a torque wrench to seat Locator Abutments. 15mm in length.


LSDT-21MM **Locator Square Drive Tool (21mm)**
Used with a torque wrench to seat Locator Abutments. 21mm in length.



products shown not to scale

BALL ABUTMENTS AND COMPONENTS

Ball Abutments

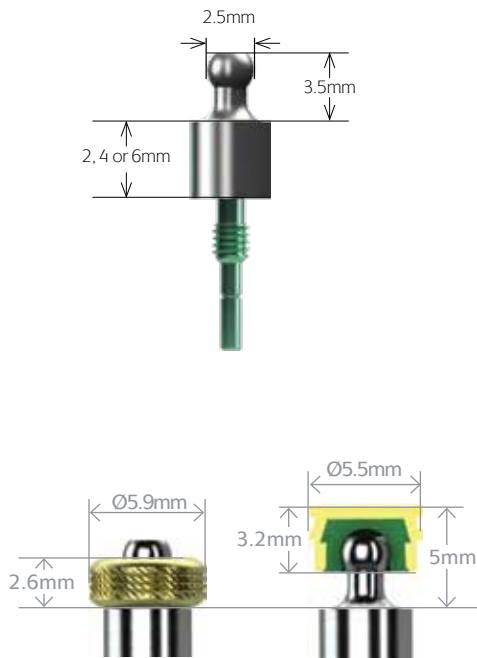

263-002 3.5mm Ball Abutment, 2mm Collar

263-004 3.5mm Ball Abutment, 4mm Collar

263-006 3.5mm Ball Abutment, 6mm Collar

264-002 4.0mm Ball Abutment, 2mm Collar

264-004 4.0mm Ball Abutment, 4mm Collar

264-006 4.0mm Ball Abutment, 6mm Collar

265-002 5.0mm Ball Abutment, 2mm Collar

265-004 5.0mm Ball Abutment, 4mm Collar

265-006 5.0mm Ball Abutment, 6mm Collar

Ball Abutments are used for retention of tissue-supported overdentures. They are recommended for relatively parallel implants. See below and on the facing page for attachment options. Ball Abutments may be used for indirect transfer impressions. Ball Abutment Analogs on facing page are used for producing a working cast. Seat with .050" (1.25mm) Hex Driver. Final torque: 30 Ncm. Titanium Alloy.

O-ring Attachment Set

Standard O-ring attachment for processing into denture.



260-100 O-ring Attachment Set

Includes: (1) O-ring encapsulator, (2) Processing O-rings and (2) Clinical O-rings.

O-ring Individual Components



260-300 O-ring Encapsulator
Female receptacle processed into denture. Titanium. 2 per package.



260-220 Processing O-ring
Used for lab processing applications. Buna. 12 per package.



260-210 Clinical O-ring
Used for clinical applications. Silicone. 12 per package.

BALL ABUTMENT COMPONENTS

Ball Attachment Set

BCAS Ball Attachment Set

Includes: (1) Titanium Housing, (3) Female Nylon Inserts – white (more retention), pink (less retention), black (lab processing) and (1) Protective Disk (ref. BCPD, protects tissue during impression making or denture pick-up)

The Ball Attachment system offers several advantages over traditional O-ring attachments. It requires less mesial/distal/buccal/lingual space and offers four different levels of retention. The nylon retention inserts can be changed chairside.



Ball Abutment Components

BCAHT Attachment Housings - Titanium

For Resin pickup or Soldering.
2 per package.



BCIY Yellow Nylon Insert

Clinical use. 2 per package.
Extra Soft Retention: 500-550g



BCIB Black Nylon Insert

Lab Processing and Chair-side
Denture Pick-up. 2 per package.



BCIP Pink Nylon Insert

Clinical use. 2 per package.
Soft Retention: 800-950g



BCIW White Nylon Insert

Clinical use. 2 per package.
Standard Retention: 1200-1300g



BCIG Green Nylon Insert

Clinical use. 2 per package.
Elastic Retention.



BCDR Directional Rings

Used for obtaining parallelism.
0°, 7° and 14° rings. Set of 3.



BCIST Insert Seating Tool

Used to seat nylon inserts in attachment housings.



BCR Reamer

Used to adjust retention of nylon inserts.



Ball Abutment Analogs



260-400

3.5/4.0mm Ball Abutment Analog (2 pack)


260-500

5.0mm Ball Abutment Analog (2 pack)

Used to represent the Ball Abutment/Implant assembly in the working cast. Only used in conjunction with Ball Abutments. Titanium Alloy.

ABUTMENT FOR SCREW COMPONENTS


Abutment for Screw - Non-hexed

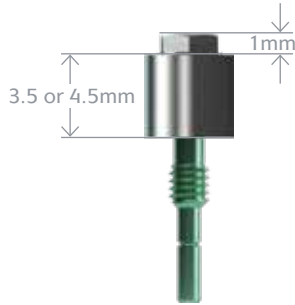

AFS353	3.5mm Abutment for Screw, 3.5mm collar
AFS354	3.5mm Abutment for Screw, 4.5mm collar
AFS403	4.0mm Abutment for Screw, 3.5mm collar
AFS504	4.0mm Abutment for Screw, 4.5mm collar
AFS503	5.0mm Abutment for Screw, 3.5mm collar
AFS504	5.0mm Abutment for Screw, 4.5mm collar

Used for multiple-unit restorations including: screw-retained restorations at the abutment level, cast alloy bars for overdentures and fixed/detachable (hybrid) restorations. Comes packaged with the Cover Cap (PXABCC). Titanium Alloy.

Angled Abutment , Abutment for Screw - Hexed

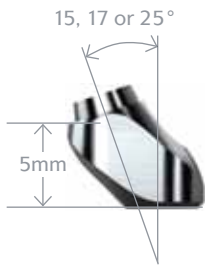

AFS355-17	3.5mm Angled Abutment for Screw, 17 Degree
AFS405-15	4.0mm Angled Abutment for Screw, 15 Degree
AFS405-25	4.0mm Angled Abutment for Screw, 25 Degree
AFS505-15	5.0mm Angled Abutment for Screw, 15 Degree

Comes packaged with the Cover Cap (PXABCC) and Abutment Screw (130-300 or 140-300). Final torque: 30Ncm. Titanium Alloy.


Gold & Plastic Copings, Abutment for Screw - Non-hexed



PXABGC

Gold Coping, Abutment for Screw

May be trimmed for height. Packaged with one Screw (regular), Abutment for Screw (see PXABS). Gold Alloy base with natural acetyl (Delrin® or Pomalux®).

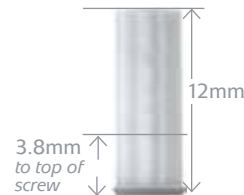

PXABPC

Plastic Coping, Abutment for Screw

May be trimmed for height. Packaged with one Screw (regular), Abutment for Screw (see PXABS). Natural acetyl (Delrin® or Pomalux®).

Instant Fixed Overdenture Abutment - Non-hexed



PXIFO

Instant Fixed Overdenture Abutment

Used with Abutments for Screw for fabrication of immediate provisional overdentures. See Fixed Provisional Overdenture (L0155) for more details. Packaged with one Screw (regular), Abutment for Screw (see PXABS). Abutment height is 9mm. Final torque: 30Ncm. Titanium Alloy.

Cover Cap, Abutment for Screw



PXABCC

Cover Cap, Abutment for Screw

Packaged with all Abutments for Screw. Hand-tighten with .050" (1.25mm) Hex Driver. Titanium Alloy.

*Abutment for Screw components are now universal for Internal, Tapered Internal, External and Single-stage provided the new Abutments for Screw are used.
For legacy External Abutment for Screw component availability, contact Customer Care.

products shown not to scale


ABUTMENT FOR SCREW COMPONENTS

Direct Pick-up Coping, Abutment for Screw

PXABDC

Direct Pick-up Coping, Abutment for Screw

Used to make a direct pick-up impression (open-tray) at the abutment level. Packaged with the Screw (long), Abutment for Screw (PXABSL). Used only with Abutment for Screw. Titanium Alloy.

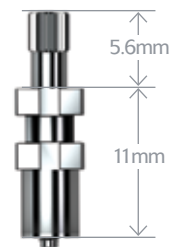


Indirect Transfer Coping, Abutment for Screw

PXABIC

Indirect Transfer Coping, Abutment for Screw

Used to make an indirect transfer (closed-tray) impression at the abutment level. Used only with Abutment for Screw. Titanium Alloy.

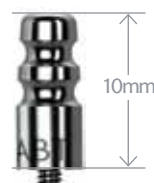

Analog, Abutment for Screw

PXABA

Analog, Abutment for Screw

Used to represent the Abutment for Screw/Implant assembly in the working cast. Not for use with implant-level impressions.



Screw (regular or long), Abutment for Screw

PXABS

Screw (regular), Abutment for Screw (5 pack)

Used to retain bars or prostheses fabricated with the Abutment for Screw Copings. Utilizes the .050" (1.25mm) Hex Driver. Final torque: 30 Ncm. Titanium Alloy.

PXABSL

Screw (long), Abutment for Screw (5 pack)

Used in the lab when a longer screw is desired. Used only with the Abutment for Screw. Up to 7mm can be removed from the screw without losing the hex engagement. Hand-tighten or torque to 30Ncm depending on application. Titanium Alloy.



Hex Adapter, Abutment for Screw

PXHA

Hex Adapter, Abutment for Screw*

Required for placement of Abutment for Screw. May be driven by either Hand Wrench, Torque Wrench or AS123 Hand Unit.



*instrument o-rings & c-rings wear out over time. If an instrument is no longer held securely by its associated driver, order a replacement ring through Customer Care.

products shown not to scale

Prosthetic Instrumentation System



320-000

Complete Prosthetic Instrumentation System (shown)

Includes: AS123 Hand Unit / Hand Wrench / 30 Ncm Torque Wrench / .050" (1.25mm) Hex Driver, Regular / .050" (1.25mm) Hex Driver, Long / Hex Adapters for Abutment for Screw / 4mm Square Drive Extender / Sterilization Tray.

320-101

Basic Prosthetic Instrumentation Kit (not shown)

Identical to the Complete Prosthetic Instrumentation System (see above), but excludes the AS123 Hand Unit and Hex Adapters for Abutment for Screw.

300-070

Prosthetic Sterilization Tray (not shown)

Autoclavable tray for prosthetic instrumentation (included with the 320-000 and 320-101).



300-400

Hand Wrench*

Used on drive end of AS123 Hand Unit. Also fits individual Hex Drivers / Adapters and Bone Taps.



300-206

4mm Square Drive Extender*

Replaced 300-205 starting in June 2010. Includes PEEK C-ring for durable retention in Ratchet. Cannot be used with bone taps.



300-100

AS123 Hand Unit*

Provides improved vision and easy access to prosthetic components in posterior regions of the mouth. Hand Wrench and Drivers must be purchased separately. May also be purchased as part of the 320-000 Complete Prosthetic Instrumentation System (see above).



300-430

30 Ncm Torque Wrench

Used to deliver torque to prosthetic components.



300-350 (regular)

300-351 (long)

.050" (1.25mm) Hex Driver Regular or Long

Used to drive all surgical and prosthetic components with the exception of the Abutment for Screw.



PXHA

Hex Adapter, Abutment for Screw*

Required for placement of Abutment for Screw. May be driven by either Hand Wrench, Torque Wrench or AS123 Hand Unit.

*instrument o-rings & c-rings wear out over time. If an instrument is no longer held securely by its associated driver, order a replacement ring through Customer Care.
products shown not to scale

ANCILLARY PROSTHETIC INSTRUMENTATION

IMPAH

Abutment Clamp

Used to hold two-piece abutments during delivery and tightening of the Abutment Screw.



Abutment Clamp

ATW

ITL Precise Adjustable Torque Wrench

Designed to place both implants and abutments with 9 distinct torque settings (15, 20, 25, 30, 35, 40, 45, 50 and 60 Ncm). A simple twist of the handle locks in precision-engineered torque values and guarantees accuracy and repeatability.

Adjustable Torque Wrenches



AGYR-15500

Anthogyr Torq Control® Torque Wrench

Ergonomic design is the ideal solution for access to screws placed in the posterior. The 7 predetermined torque values (10, 15, 20, 25, 30, 32 and 35 Ncm) make it a tool of extreme precision.



293-100

3.5mm Implant Analog Handle

294-100

4.0mm Implant Analog Handle

295-100

5.0/6.0mm Implant Analog Handle

Used to hold External abutments for chairside or laboratory preparation, these handles mimic the implant/analog hex geometry. Abutments are secured to the handle with the standard Abutment Screw. Comes in three platform sizes: 3.5, 4.0mm and 5.0/6.0mm.

Abutment Prepping Handles



135-351

.050" (1.25mm) One-piece Hex Driver

134-350

.050" (1.25mm) Handpiece Hex Driver

134-450

.050" (1.25mm) Handpiece Hex Driver, Long

For installation and removal of Cover Screws, Healing Abutments and Abutment Screws. The Handpiece Hex Drivers are used with latch-type contra-angle handpieces. The Handpiece Hex Driver, Long (134-450) is 5mm longer than the standard version (134-350).

Hex Drivers



123-001

3.5mm Implant Clean-out Tap Tool

124-001

4.0/5.0/6.0mm Implant Clean-out Tap Tool

122-170

Abutment for Screw Clean-out Tap Tool

Used to re-thread External implants or Abutments for Screw where the internal threads have become damaged. Requires a standard surgical Ratchet (130-000) or Hand Wrench (300-400) as a drive mechanism.

Clean-out Tap Tools



**instrument o-rings & c-rings wear out over time. If an instrument is no longer held securely by its associated driver, order a replacement ring through Customer Care.
products shown not to scale*



This Surgical Manual serves as a reference for using the External implants and surgical instruments. It is intended solely to provide instructions on the use of BioHorizons products. It is not intended to describe the methods or procedures for diagnosis, treatment planning, or placement of implants, nor does it replace clinical training or a clinician's best judgment regarding the needs of each patient. BioHorizons strongly recommends appropriate training as a prerequisite for the placement of implants and associated treatment.

The procedures illustrated and described within this manual reflect idealized patient presentations with adequate bone and soft tissue to accommodate implant placement. No attempt has been made to cover the wide range of actual patient conditions that may adversely affect surgical and prosthetic outcomes. **Clinician judgment as related to any specific case must always supersede any recommendations made in this or any BioHorizons literature.**



Before beginning any implant surgical procedure with BioHorizons implants:

- Read and understand the Instructions for Use that accompany the products.
- Clean and sterilize the surgical tray and instruments per Instructions for Use.
- Become thoroughly familiar with all instruments and their uses.
- Study Surgical Kit layout and iconography.
- Design a surgical treatment plan to satisfy the prosthetic requirements of the case.

Indications for Use

BioHorizons dental implants may be used in the mandible and maxilla for use as an artificial root structure for single tooth replacement or as abutments for bridgework and denture retention.

VIP Treatment Planning



Virtual Implant Placement (VIP) treatment planning software is a user-friendly solution that reduces clinical challenges and enhances post-operative outcomes.

- Interactive 2D and 3D treatment planning
- Self processing DICOM converter
- Case viewer available for download from BioHorizons website

Two-stage Protocol

Two-stage surgery was the original protocol developed for placing modern dental implants. In this type of surgery, the implant is placed below the soft tissue and protected from occlusal function and other forces during osseointegration. A low-profile Cover Screw is placed on the implant to protect it from the ingress of soft tissue.

Following osseointegration, a second surgery exposes the implant and a transmucosal Healing Abutment is placed to allow for soft tissue healing and development of a sulcus. Prosthetic restoration begins after soft tissue healing.



External implant with Cover Screw in a two-stage protocol.

Single-stage Protocol

Single-stage surgery may be accomplished with the External implant by placing a healing abutment at the time of initial implant surgery. This eliminates the need for a second surgery to expose the implant. Although the implant is not in occlusal function, some forces can be transmitted to it through the exposed transmucosal element.

Prosthetic restoration begins following osseointegration of the implant and soft tissue healing.



External implant with removable Healing Abutment in a single-stage protocol.

Non-functional Immediate Restoration

Single-stage surgery with non-functional immediate provisionalization provides the patient a non-functioning provisional prosthesis early in the treatment plan. An abutment is placed on the implant at or shortly after surgery, and a provisional restoration is secured to it with temporary cement. The provisional can help contour the soft tissue profile during healing.



External implant restored with a non-functional provisional prosthesis.

Immediate Function Restoration

Single-stage surgery with immediate function is possible in good quality bone where multiple implants exhibiting excellent initial stability can be splinted together. Splinting implants together can offer a significant biomechanical advantage over individual, unsplinted prostheses.



External implants with a splinted prosthesis in immediate function.


products shown not to scale



During implant placement, clinicians must apply their best judgment as to the appropriate spacing for individual patient conditions.

Spacing considerations for BioHorizons External implants (measurements are taken at the osseous crest):


The osteotomy centerpoint required to maintain a 1.5mm implant-to-tooth spacing (generally accepted) is $\frac{1}{2}$ [implant body diameter] + 1.5mm.

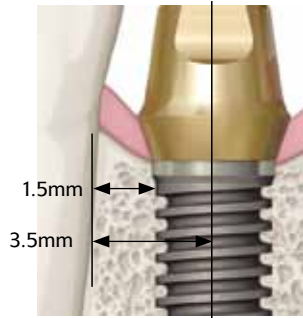

Osteotomy center 3.5mm from adjacent tooth
(4.0mm implant pictured)

body diameter	osteotomy center from adjacent tooth
3.5mm	3.3mm
4.0mm	3.5mm
5.0mm	4.0mm
5.0mm	4.5mm

The osteotomy center-to-center measurement required to maintain a 3.0mm edge-to-edge spacing (generally accepted) between External implants is $\frac{1}{2}$ [sum of 2 implant body diameters] + 3.0mm.

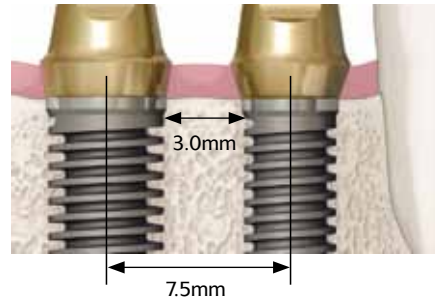

Measurement is dependent on the two implant body diameters.
(4.0 & 5.0mm implants pictured)

body diameter	3.5mm	4.0mm	5.0mm	6.0mm
3.5mm	6.5mm			
4.0mm	6.8mm	7.0mm		
5.0mm	7.3mm	7.5mm	8.0mm	
5.0mm	7.8mm	8.0mm	8.5mm	9.0mm

Implant Spacer / Depth Probe

Purpose: Multi-function instrument for intraoral measurements.


- Five centimeter graduated ruler on shaft




The rectangular end of the tool provides intraoral measurements.



Useful for marking center-to-center implant spacing on the ridge.

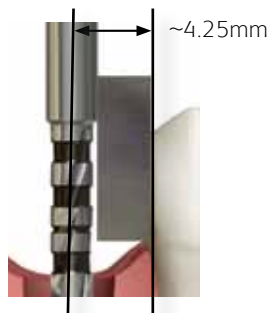

The rectangular end against an existing crown places the osteotomy ~4.25mm from the contact.

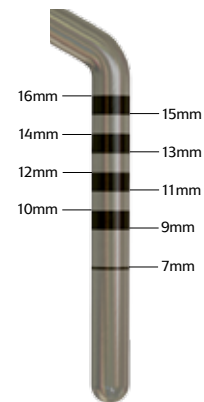


Probe tip measures osteotomy depth.
Note: these markings are different than the tapered drill markings

SURGICAL KIT LAYOUT


The External Surgical Kit uses an intuitive layout to guide the surgeon through the instrument sequence. The sequence begins in the upper left hand corner and works left-to-right and then down. Color-coded lines, instruments and grommets are matched with each implant and further aid in instrument selection and identification.

Prior to use, clean and sterilize the surgical tray and instruments per appropriate Instructions for Use and study the Surgical Kit layout, color-coding and iconography. Surgical assistants should also be thoroughly familiar with all instruments and their uses.

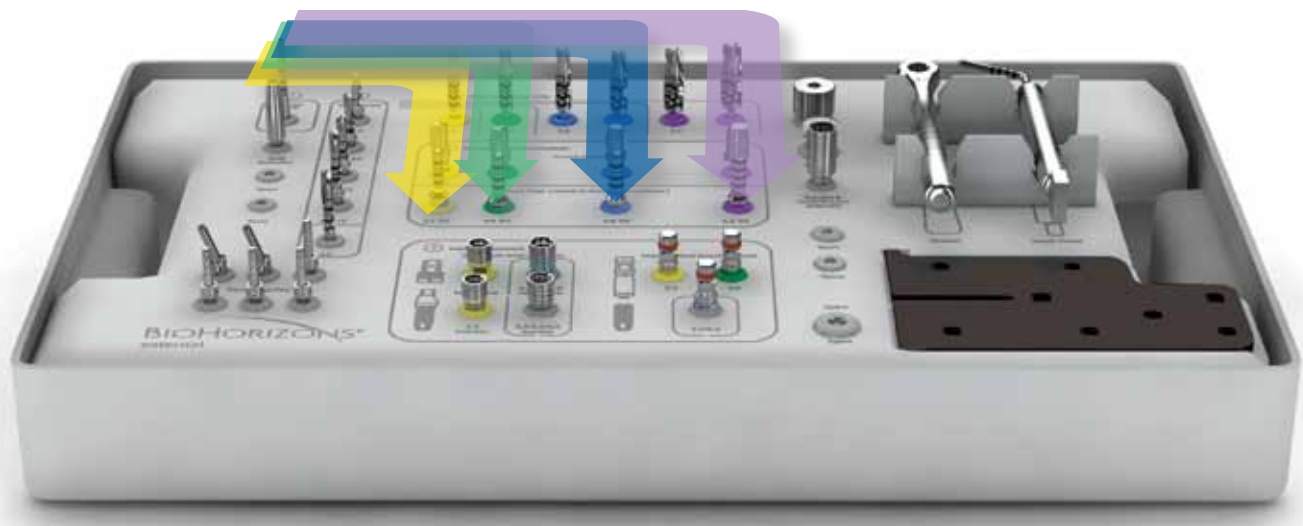

Depth Drill Sequence	Width Increasing Drill Sequence	Crestal Bone Drills	Bone Taps	Abutment-level and/or Implant-level Drivers

The Surgical Drills included in the External Surgical Kit are externally irrigated and designed to be used at speeds of 850-2,500 rpm⁹ with steady sterile irrigation. Reduced drill speed may be desired in softer bone or as drill diameter increases.

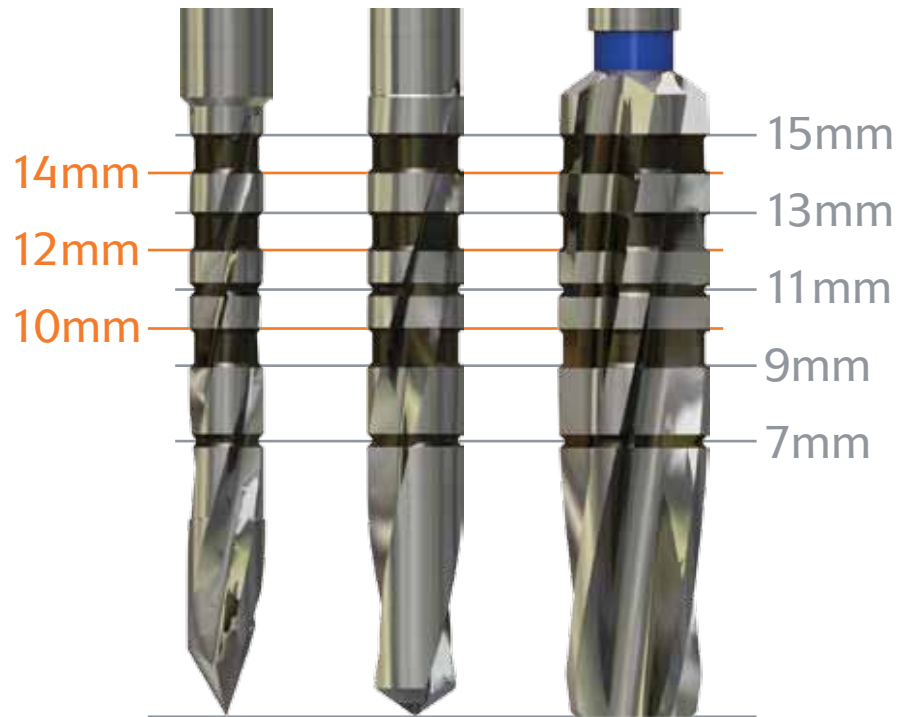

The depth marks are consistent throughout the Starter Drill, Depth Drills and Width Increasing Drills. The thick bands are each one millimeter wide. Thin lines are used to mark 7mm and 11mm.

Drilling Considerations

Peri-operative oral rinses with a 0.12% Chlorhexidine Digluconate solution have been shown to significantly lower the incidence of post-implantation infectious complications.¹⁰ A preoperative 30-second rinse is recommended, followed by twice daily rinses for two weeks following surgery.

Drilling must be done under a constant stream of sterile irrigation. A pumping motion should be employed to prevent overheating the bone. Surgical drills and taps should be replaced when they are worn, dull, corroded or in any way compromised. BioHorizons recommends the replacement of drills after 12 to 20 osteotomies.¹¹ A Drill-usage Tracking Chart is available from BioHorizons to aid offices in archiving this important information.

There is a risk of injury to the mandibular nerve associated with surgical drilling in lower posterior regions. To minimize the risk of nerve injury, it is imperative the clinician understand the drill depth markings in order to correlate implant length with the actual drilling depth to produce the desired vertical placement of the implant.

OSTEOTOMY INITIALIZATION

2.0mm Starter Drill

Purpose: Initiates osteotomy.


- Chisel-tip design eliminates “skating” on osseous crest
- Initiates osteotomy to desired depth
- Prepares site for Paralleling Pins



2.5mm Depth Drill

Purpose: Sets osteotomy depth following use of the 2.0mm Starter Drill.


- Efficient cutting drill design collects bone for autografting

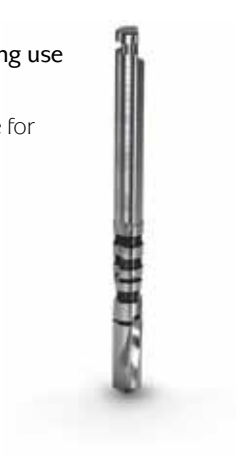

2.5mm Depth Drill with Stop

Purpose: Sets osteotomy depth following use of the 2.0mm Starter Drill.

- Fixed circular ring acts as a definitive drill stop
- One drill length for each implant length
- 1 mm laser-etched line guides supracrestal implant placement
- BioHorizons Surgical Kit includes spare slots for Depth Drills with Stops
- Please call for availability

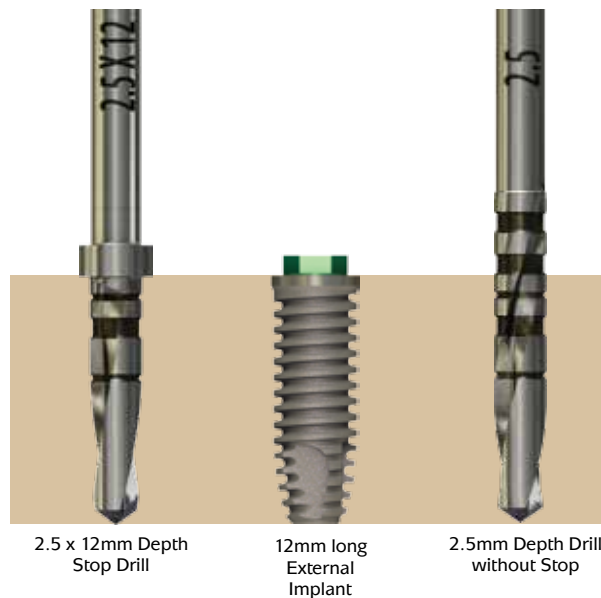

2.5 x 12mm Depth
Stop Drill

12mm long
External
Implant

2.5mm Depth Drill
without Stop

products shown not to scale

Parallel Pins

Purpose: Evaluation of osteotomy position and angle.


- Provided straight and with a 20° angle
- Used after 2.0mm Starter Drill and 2.5mm Depth Drills
- 9mm shank for radiographic evaluation of proximity to adjacent anatomy
- Hub diameter is 4.0mm



Width Increasing Drills

Purpose: Incrementally widens the osteotomy without creating excessive heat.

- Depth-marked for reference. Non-end cutting geometry for added safety
- Efficient cutting geometry collects bone for autografting
- Color-coded by the smallest diameter implant body with which it is used

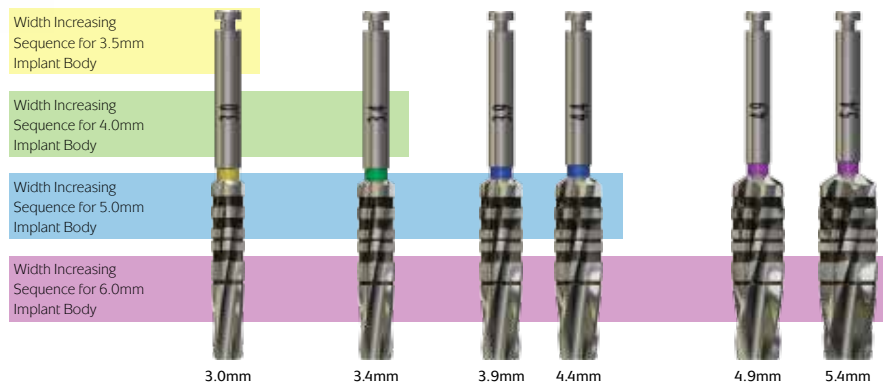

Variations in bone density or possible voids require surgeons to observe the drill's depth markings as the primary determinant of drilling depth, rather than relying exclusively on the non-end cutting geometry to stop the drill.



If a Width Increasing Drill fails to reach the planned depth as shown by the depth marks, use the 2.5mm Depth Drill to re-establish the depth and recommence the Width Increasing Sequence with the 3.0mm Width Increasing Drill.


CRESTAL BONE DRILLS & BONE TAPS

Crestal Bone Drills

Purpose: Removes cortical bone at ridge crest to facilitate pressure-free seating of the implant collar.

- Site Specific. Indicated when dense cortical bone is present at crest
- Rounded non-end cutting hub centers drill in osteotomy
- Used following the final Width Increasing Drill for each implant
- 850-2,500 rpm with steady irrigation¹¹

Do not fully seat the drill if all or a portion of the collar is to be left supracrestal. Fully seating the drill allows the implant to be placed with the prosthetic platform level with the crestal ridge.

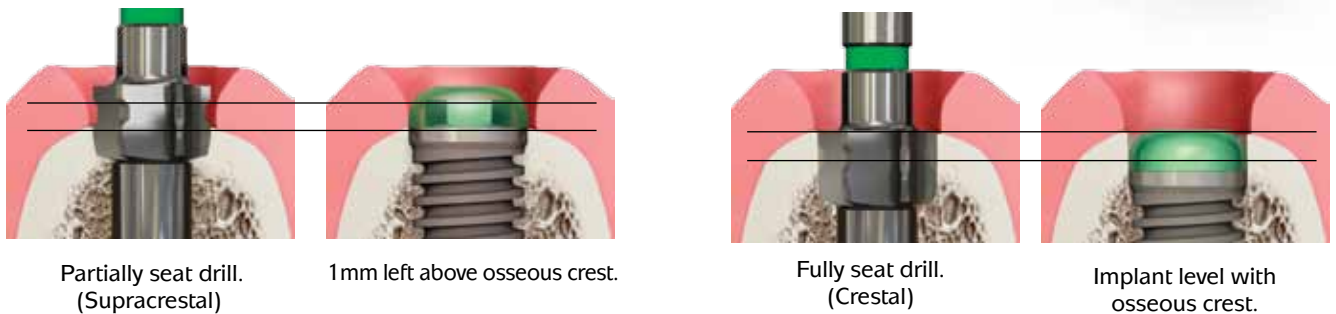

- Preparing to Level 1 allows placement of the prosthetic platform level with the crest.
- Preparing to Level 2 allows placement of the Cover Screw level with the crest, *provided the osteotomy has been prepared to accommodate the additional depth.*

Bone Taps

Purpose: Prepares dense cortical bone for implant threads.

- Site specific
- 30 rpm or less¹²
- Final instrument prior to implant placement
- Square drive shaft interfaces with Ratchet and Hand Wrench



Place the tip of the Bone Tap into the osteotomy, apply firm apical pressure and begin rotating slowly in a clockwise direction (30 rpm or less is recommended)¹². When the threads engage, allow the tap to feed without excessive pressure. To remove, rotate the Bone Tap in a counter-clockwise direction, allowing it to back out of the osteotomy.

Do not pull on the Bone Tap to remove it from the site.



products shown not to scale

Implant Packaging

External implants are provided in a blister pack with multiple peel-and-stick labels for affixing to the patient's chart. When the lid of the blister tray is removed, the sterile implant vial is exposed and may then be placed in the sterile field. While holding the vial in an upright fashion, remove the cap by rotating it in a counter-clockwise direction. The implant can then be removed from the vial by engaging the premounted *3inOne* Abutment with the appropriate Adapter.



Handpiece and Ratchet Adapters (Abutment-level)

Purpose: Engages the *3inOne* Abutment allowing the implant to be driven into the osteotomy.


- Premounted *3inOne* Abutment serves as the surgical drive mount
- O-ring secures implant while it is carried to osteotomy
- Electric handpiece or manual insertion options
- 30 rpm or less¹²



Handpiece and Ratchet Adapters engage External implants via the square on the coronal aspect of the *3inOne* abutment. Remove the cap from the inner vial and seat the chosen Adapter. Remove the implant from the vial and carry it to the osteotomy on the Adapter, taking precautions not to touch the implant surface during the transfer.

Place the apex of the implant into the osteotomy, apply firm apical pressure and begin rotating in a clockwise direction (30 rpm or less). Once the apical threads engage the bone, allow the implant to feed without excessive pressure.

Overtightening the implant in the osteotomy may cause osseous microfracture or pressure-induced necrosis of the adjacent bone. Manual seating via the Adapter for Ratchet may be desired to gain a tactile sense of final implant placement. If too much resistance is felt during insertion, remove the implant and revise the osteotomy with the appropriate Crestal Bone Drill or Bone Tap as deemed necessary to reduce insertion torque.

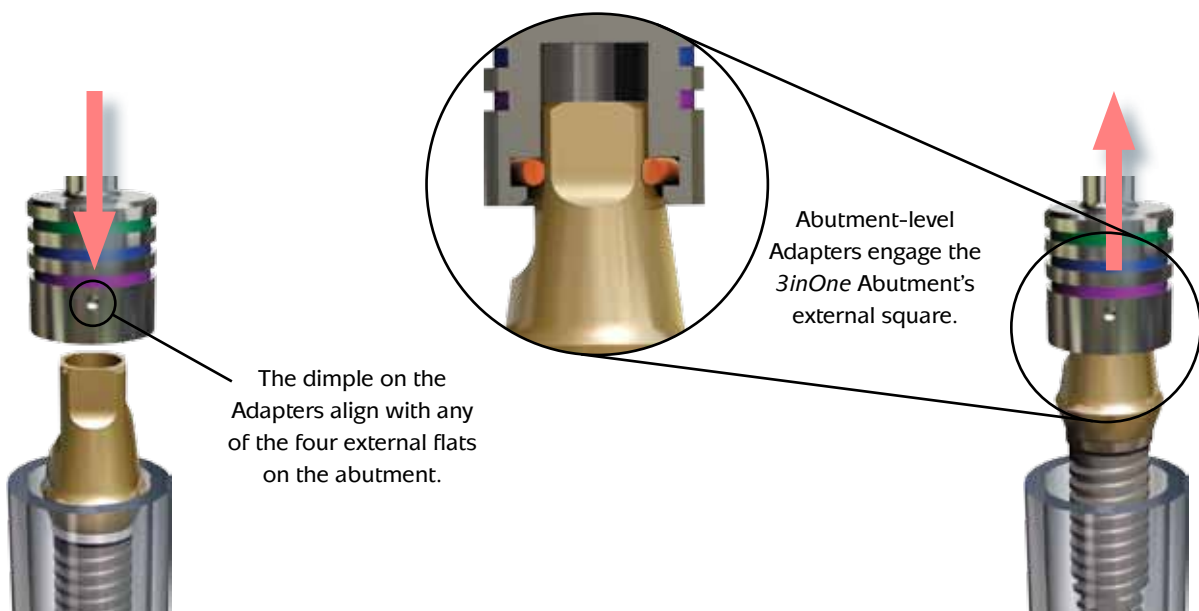

3inOne Abutment Removal

To remove the *3inOne* Abutment, engage the Abutment Screw with the .050" (1.25mm) Hex Driver. Apply firm apical pressure to the Hex Driver and rotate counter-clockwise until the screw is completely disengaged from the implant body.

In soft bone, or when the implant lacks initial stability, an Abutment Clamp (ref. IMPAH, sold separately) should be used to grasp the outside of the abutment to provide counter-torque during the loosening of the Abutment Screw.

The *3inOne* Abutment and the Abutment Screw should be retained with the patient's chart. They can later be used in the impression making procedure and as a temporary or final abutment for cement retention.



Hex Orientation

The longest flat surface on the external aspect of the *3inOne* Abutment is indexed to one of the six flats of the implant's external hexagon. In most cases one of the hex flats should be oriented to the facial aspect, as it allows for angulation correction with stock angled abutments. Placing the long flat of the *3inOne* Abutment to the facial also leaves more room for porcelain in that area on the final prosthesis.



The implant's rotational position can be adjusted following removal of the *3inOne* Abutment using the Implant-level Insertion Tools.



Insertion Tools (Implant-level)

Purpose: Engages external hex allowing implant's position to be adjusted in the osteotomy.




- May be used following removal of the *3inOne* Abutment
- Offers a narrower path of insertion and better clearance than the *3inOne* Abutment / Abutment-level Adapter option
- Square drive shaft interfaces with the Ratchet and Hand Wrench

Mate the appropriate Insertion Tool with the external hex of the implant and lightly screw into place. Engage the Insertion Tool with either the Ratchet or Hand Wrench and rotate the implant to the required position.



COVER SCREWS & HEALING ABUTMENTS

Cover Screws

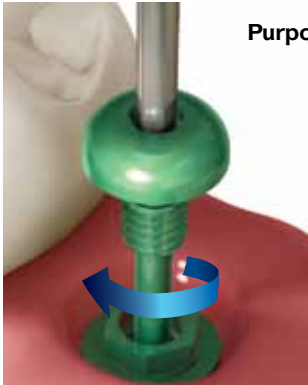

Purpose: Protects prosthetic platform in two-stage (submerged) surgical protocols.

- Packaged with each implant in the sterile inner vial
- Hand-tighten (10-15 Ncm)¹³ utilizing .050" (1.25mm) Hex Driver
- Color-coded by implant body diameter



Remove the *3inOne* Abutment and thoroughly irrigate the inside of the implant to remove blood and other debris. Unscrew the Cover Screw (included in each implant vial) from its holder and screw it into the implant using the .050" (1.25mm) Hex Driver. Following placement of the Cover Screw, the surgical site should be irrigated and the soft tissue adapted in a normal surgical fashion. Take precautions to prevent the Cover Screw from being aspirated by the patient.

Healing Abutments

Purpose: Transmucosal element for single-stage surgical protocol or for soft tissue healing period following second-stage uncover.



- Select by height and emergence profile
- Hand-tighten (10-15 Ncm)¹³ utilizing .050" (1.25mm) Hex Driver
- Color-coded by implant body diameter
- Encoded for easy intraoral identification, for example:
4 x 4.5 indicates 4.0mm implant diameter x 4.5mm high



Healing Abutments are placed after uncover in a two-stage surgical protocol, or in lieu of a Cover Screw in a single-stage (non-submerged) protocol. Prior to seating the Healing Abutment, thoroughly irrigate the inside of the implant to remove blood and other debris. An antibacterial paste may be placed on the screw portion to decrease the risk of bacterial growth within the implant body during the healing phase. Following seating, irrigate the surgical site and adapt the soft tissue in normal surgical fashion. A gingivectomy or apically positioned flap technique may be used to reduce the soft tissue thickness and to decrease sulcular depth around the implant. The suture groove on the Healing Abutment may be used to apically position the soft tissue flap. Take precautions to prevent the Healing Abutment from being aspirated by the patient.



All 7mm implants come packaged with a special Cover Screw and Abutment Screw that have shorter guide pins for complete seating into the 7mm implant. Standard length Abutment Screws and Healing Abutments can be used provided the guide pin is carefully removed at the notch using a high-speed carbide disk. Please refer to the External -7mm Length Implants manual (ref. L0151) for detailed information.

Bone Profiler



Purpose: Remove and contour excess bone and soft tissue from the area of the prosthetic platform.

- Compatible with latch-type, speed-reducing handpieces
- 850-2,500 rpm drill speed with steady sterile irrigation⁹
- Cover Screw protects implant platform
- Bone Profiler cuts away excess bone and soft tissue
- Color-coded by specific prosthetic platform



Do not use the Profiler without the Cover Screw in place.



The Cover Screw aligns the Bone Profiler and protects the implant from damage. The Profiler is used in a latch-type, reduction handpiece under copious amounts of sterile irrigation. Following removal of the excess bone and soft tissue, unscrew the Cover Screw from the implant and seat the desired prosthetic element.

Bur Testing Gauge

Also called a “Go / No-Go Gauge,” the Bur Testing Gauge is used to verify the dimensional accuracy of drill shanks of latch-type burs. Burs in proper condition **WILL** fit into the larger diameter hole, but **WILL NOT** fit into the smaller hole (marked red).

Burs that fail either of these criteria are unfit for use and may become stuck in the handpiece if used. The gauge is included with all W&H starter packages, and may also be ordered separately.



Surgical Kit Cleaning

All BioHorizons Surgical Kits are provided non-sterile and must be cleaned and sterilized prior to use following the associated Instructions for Use. Always remove instruments from packaging prior to sterilization, and remove and discard packaging materials used to stabilize and secure kits during shipment. Double-check all surgical instruments to ensure their functionality prior to surgery. Backup sterile drills are also recommended.



Caution: The use of hydrogen peroxide or other oxidizing agents will cause damage to the surface of the instruments. Towel- or air-dry all instrumentation before sterilizing. Drills should be replaced after approximately 12 to 20 osteotomy cycles, depending on the bone density.¹¹

Proper testing, cleaning and calibration of sterilization equipment should occur frequently to assure that the units are in proper working order. Equipment operating conditions vary and it is the responsibility of each dental office to ensure that the proper sterilization technique for instrumentation is followed.

Post-operative Considerations

A period of unloaded healing is often recommended, depending on individual patient healing rates and initial implant stability. Each case must be independently evaluated. This unloaded healing period allows for integration between the bone and implant surface.

The patient must be instructed to follow a post-surgical regimen including cold packs for 24 hours post-implantation. The patient's diet should consist of soft foods and possibly liquid dietary supplements. Pharmacological therapy should be considered as the patient's condition dictates.

If a removable prosthesis is used during the initial healing phase, it is recommended that a soft liner material be used to prevent pressure on the surgical site. This soft liner should be relieved over the implant site. The patient should be checked periodically to monitor healing of the soft tissues and bone using clinical and radiographic evaluations.

Ongoing hygiene for the implant patient is vital. Hygiene recall appointments at three month intervals are suggested. Instruments designed for implant scaling, such as Implacare® instruments from Hu-Friedy® should be utilized. The stainless steel handles may be fitted with assorted tip designs used for hygiene on natural teeth. The Implacare® scalers will not damage implant abutments and contain no glass or graphite fillers that can scratch titanium implant abutments.

Symbol Descriptions for Product Labeling

The diagram shows a central image of a BioHorizons product label for a 4.0 x 12 External Implant. Callouts point to various elements on the label and provide their meanings:

- REF**: Reference/article number
- LOT**: Lot/batch number
- Hourglass icon**: Use before expiration date (YYYY-MM)
- Calendar icon**: Manufacture date (YYYY-MM)
- STERILE IR**: Sterile by gamma irradiation
- NON-STERILE**: Non-sterile
- Rx Only**: Caution: Federal (USA) law restricts these devices to the sale, distribution and use by, or on the order of, a dentist or physician.
- Artwork label number**: 0473
- Prosthetic platform**:
 - 3.5mm Prosthetic Platform (Yellow)
 - 4.0mm Prosthetic Platform (Green)
 - 5.0mm Prosthetic Platform (Blue)
 - 6.0mm Prosthetic Platform (Purple)
- Single use only**: Indicated by a crossed-out circular arrow icon.
- Refer to Instructions for Use**: Indicated by a triangle with an exclamation mark icon.
- CE 0473**: BioHorizons products carry the CE mark and fulfill the requirements of the Medical Devices Directive 93/42/EEC.
- EU Authorised Representative**: Quality First International, London E7 0QY United Kingdom. Tel. & Fax: +44-208-522-1937.

The label itself contains the following text: 4.0 x 12 External Implant, BIOHORIZONS Birmingham, AL 35244 USA, REF 4012D3, RBT, D3 Threadform 4.0 x 12mm, LOT YYXXXXX, YYY-MM expires, YYY-MM manufacture date, STERILE IR gamma irradiated, Rx Only, L4012D3 Rev A, BioHorizons External REF 4012D3 LOT YYXXXXX 4.0 x 12mm, BioHorizons External REF 4012D3 LOT YYXXXXX 4.0 x 12mm, and the CE mark with number 0473.

References

1. Please see BioHorizons document #ML0130.
2. Functional Surface Area: Thread-Form Parameter Optimization for Implant Body Design. J. Todd Strong, Carl E. Misch, DDS, MDS, Martha W. Bidez, PhD, Prasad Nalluri. *Compendium / Special Issue*. Vol 19, No. 3, 1998.
3. Histologic and Histomorphometric Findings from Retrieved, Immediately Occlusally Loaded Implants in Humans. Romanos GE, Testori T, Degidi M, Piattelli A. *J Periodontol*. 2005 Nov;76(11):1823-1832.
4. Effects of Implant Thread Geometry on Percentage of Osseointegration and Resistance to Reverse Torque in the Tibia of Rabbits. Steigenga J, Al-Shammari K, Misch C, Nociti Jr. FH, Wang H-L. *J Periodontol*. 2004;75(9):1233-1241.
5. A Prospective Multicenter Clinical Investigation of a Bone-Quality Based Implant System. Kline R, Hoar J, Beck G, Hazen R, Resnik R and Crawford E. *Implant Dentistry*, 224-234, Vol. 11 No. 3, 2002
6. Five-Year Prospective Study of Immediate / Early Loading of Fixed Protheses in Completely Edentulous Jaws with a Bone Quality-Based Implant System. Misch CE and Degidi M. *Clinical Implant Dentistry and Related Research*, 17-27, Vol. 5 No. 1, 2003
7. Immediate Functional and Non-Functional Loading of Dental Implants: a 2- to 60-Month Follow-Up Study of 646 Titanium Implants. Degidi M and Piattelli A. *J Periodontology* 2003;74:225-241
8. Non-Functional Immediate Teeth in Partially Edentulous Patients: A Pilot Study of 10 Consecutive Cases Using the Maestro Dental Implant System. Misch CE. *Compendium Special Issue*, 25-36, 1998
9. Density of Bone: Effect on Surgical Approach and Healing. CE Misch. *Contemporary Implant Dentistry. Second Edition*. Mosby: St. Louis, 1999. 371-384.
10. The Influence of 0.12% Chlorhexidine Digluconate Rinses on the Incidence of Infectious Complications and Implant Success. Lambert P, Morris H and Shigeru O. *J Oral Maxillofac Surg*, 1997, 55:25-30 Suppl 5
11. Heat production by 3 implant drill systems after repeated drilling and sterilization. Chacon GE, Bower DL, Larsen PE, McGlumphy EA, Beck FM. *J Oral Maxillofac Surg*. 2006 Feb;64(2):265-9
12. Root Form Surgery in the Edentulous Mandible: Stage I Implant Insertion. CE Misch. *Contemporary Implant Dentistry. Second Edition*. Mosby: St. Louis, 1999. 347-369.
13. Implant abutment screw torque generated by general dentists using a hand driver in a limited access space simulating the mouth. E. Hill, S. Phillips and L. Breeding. *J. of Oral Implantology*. 2007 Vol. 33(5):277-279

ORDERING & WARRANTY INFORMATION

Product Support Specialist: _____

Cell phone: _____

Fax: _____

BioHorizons Lifetime Warranty on Implants and Prosthetics: All BioHorizons implants and prosthetic components include a Lifetime Warranty. BioHorizons implant or prosthetic components will be replaced if removal of that product is due to failure (excluding normal wear to overdenture attachments).

Additional Warranties: BioHorizons warranties instruments, surgical drills, taps, torque wrenches and Virtual Implant Placement (VIP) treatment planning software.

(1) Surgical Drills and Taps: Surgical drills and taps include a warranty period of ninety (90) days from the date of initial invoice. Surgical instruments should be replaced when they become worn, dull, corroded or in any way compromised. Surgical drills should be replaced after 12 to 20 osteotomies.¹¹

(2) Instruments: The BioHorizons manufactured instrument warranty extends for a period of one (1) year from the date of initial invoice. Instruments include drivers, sinus lift components, implant site dilators and BioHorizons tools used in the placement or restoration of BioHorizons implants.

(3) VIP treatment planning software: VIP treatment planning software warranty extends for a period of ninety (90) days from the date of initial invoice. The warranty requires that VIP be used according to the minimum system requirements.

(4) Compu-Guide surgical templates: Compu-Guide surgical templates are distributed without making any modifications to the submitted Compu-Guide Prescription Form and VIP treatment plan ("as is"). BioHorizons does not make any warranties expressed or implied as it relates to surgical templates.

Return Policy: Product returns require a Return Authorization Form, which can be acquired by contacting Customer Care. The completed Return Authorization Form should be included with the returned product. For more information, please see the reverse side of the invoice that was shipped with the product.

Disclaimer of Liability: BioHorizons products may only be used in conjunction with the associated original components and instruments according to the Instructions for Use (IFU). Use of any non-BioHorizons products in conjunction with BioHorizons products will void any warranty or any other obligation, expressed or implied.

Treatment planning and clinical application of BioHorizons products are the responsibility of each individual clinician. BioHorizons strongly recommends completion of postgraduate dental implant education and adherence to the IFU that accompany each product. BioHorizons is not responsible for incidental or consequential damages or liability relating to use of our products alone or in combination with other products other than replacement or repair under our warranties.

Compu-Guide surgical templates are ordered under the control of a Clinician. The Clinician recognizes responsibility for use. Therefore, regardless of the real or proven damages, the liability to BioHorizons is limited to the price of the product directly related to the reason for the claim.

Distributed Products: For information on the manufacturer's warranty of distributed products, please refer to their product packaging. Distributed products are subject to price change without notice.

Validity: Upon its release, this literature supersedes all previously published versions.

Availability: Not all products shown or described in this literature are available in all countries.

BioHorizons continually strives to improve its products and therefore reserves the right to improve, modify, change specifications or discontinue products at any time.

Any images depicted in this literature are not to scale, nor are all products depicted.

Direct Offices

BioHorizons USA

888-246-8338 or
205-967-7880

BioHorizons Canada

866-468-8338

BioHorizons Spain

+34 91 713 10 84

BioHorizons UK

+44 (0)1344 752560

BioHorizons Germany

+49 7661-909989-0

BioHorizons Australia


+61 2 8399 1520

BioHorizons Chile

+56 2 361 9519

Distributors

For contact information in our 80 markets, visit www.biohorizons.com

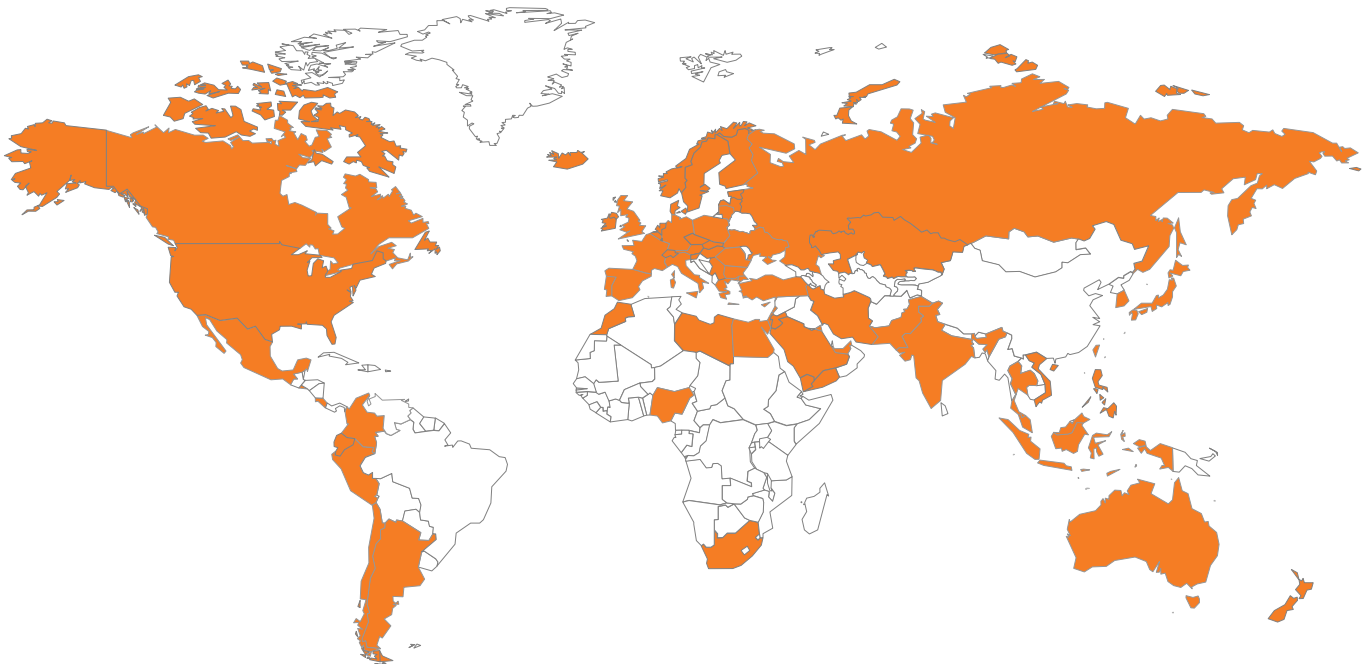

BioHorizons®, Laser-Lok®, MinerOss®, AutoTac® and Mem-Lok® are registered trademarks of BioHorizons, Inc.

Grafton® DBM and LADDEC® are registered trademarks of Osteotech, Inc. AlloDerm® and AlloDerm GBR® are registered trademarks of LifeCell™ Corporation.

Spiralock® is a registered trademark of Spiralock Corporation. Locator is a registered trademark of Zest Anchors, Inc.

Delrin® is a registered trademark of E.I. du Pont de Nemours and Company. Pomalux® is a registered trademark of Westlake Plastics Co.

Mem-Lok® is manufactured by Collagen Matrix, Inc.

Not all products shown or described in this literature are available in all countries. As applicable, BioHorizons products are cleared for sale in the European Union under the EU Medical Device Directive 93/42/EEC and the tissues and cells Directive 2004/23/EC. We are proud to be registered to ISO 13485:2003, the international quality management system standard for medical devices, which supports and maintains our product licences with Health Canada and in other markets around the globe.

Original language is English. © 2010 BioHorizons, Inc. All Rights Reserved.

shop online at
www.biohorizons.com



L0129

REV D NOV 2010